


Forensic Lab Solutions



SIRCHIE®
Forensics

www.sirchie.com

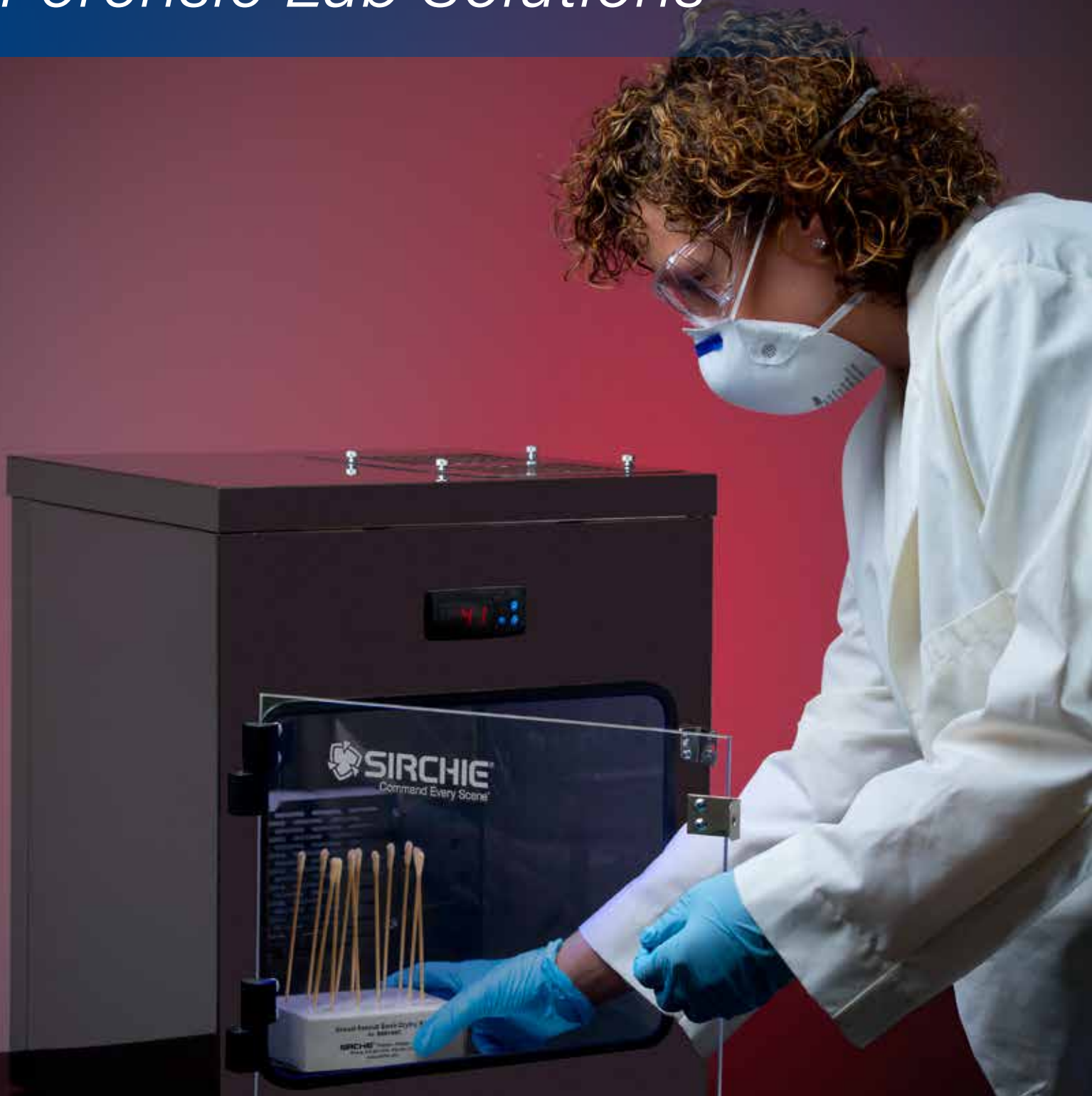

Introduction

As an investigator or a lab technician, you have a difficult and important job of finding, preserving, processing, and analyzing evidence that can solve criminal cases, bring perpetrators to justice, and ease the pain of victims' families.

We're your full-service partner in every step you take to protect, identify, document, collect, develop, examine and preserve evidence. We've worked with industry professionals to develop solutions using modern technology so you can do your job effectively and efficiently with confidence.

"Here's our Budget. Build Us a Lab."

We often hear from law enforcement professionals, crime lab directors, educators, contractors, architects, and construction companies that they need assistance with creating and building out a lab. That's where we come in. With the breadth and depth of our experience and solutions, we can guide you to the best solutions that help you meet your mission and balance your budget.



Table of Contents

Personal Protection	4-8
NARKsafe™ Protection, Packaging & Testing	9-12
Presumptive Drug Identification	13-14
Alternate Light Sources	15
Evidence Photography	16-17
Latent Fingerprint Powders	18-20
RUVIS Systems	21-25
Chemical Fingerprint Development	26-29
Ductless Drying Cabinets, Fuming Chambers & Workstations	30-33
Magnifiers & Examination Station	34-36
DNA & Biological Evidence Collection	37-42

NARKsafe™ Illicit Substance Testing Workstation

NARKsafe™ fume hood enclosures are ideal for protecting professionals during the forensic analysis of fentanyl and other dangerous substances.

See Page 12 for details!



NARKsafe NDA Narcotic Disposal Agent

NARKsafe NDA is a drug disposal solution that enables the onsite disposal of drug evidence (20 grams or less) – quickly, easily and safely, rendering the drugs unusable.

See Page 10 for details!



Forensic Swab Dryer

Use our Forensic Swab Dryer to provide a safe, controlled environment for drying swabs after specimen collection. This dryer will hold up to 72 swabs with a drying time of one hour.

See Page 41 for details!



Personal Protection

Gloves

Officials should always wear gloves when in a scene or handling an unknown substance. Ideally, you should be double-gloved in most situations. ***The DEA and the CDC recommend nitrile gloves because latex gloves are porous and could expose you to substances.***

Black Powder-Free Nitrile Gloves

These thin, comfortable, ambidextrous nitrile gloves provide extra puncture and abrasion resistance. Always use nitrile gloves while working with illicit substances and NARK presumptive drug kits.



No. SF0081S	Black, Powder-Free, Sm, 100 ea.	\$20.20
No. SF0081M	Black, Powder-Free, Md, 100 ea.	\$20.20
No. SF0081L	Black, Powder-Free, Lg, 100 ea.	\$20.20
No. SF0081XL	Black, Powder-Free, XL, 100 ea.	\$20.20

Blue Powder-Free Nitrile Gloves

These blue, powder-free, nitrile gloves provide extra puncture and abrasion resistance over latex. Sold in boxes of 100 ambidextrous gloves.



No. SF00851	Blue Powder-Free Nitrile Gloves, Sm, 100 ea.	\$20.20
No. SF00852	Blue Powder-Free Nitrile Gloves, Md, 100 ea.	\$20.20
No. SF00853	Blue Powder-Free Nitrile Gloves, Lg, 100 ea.	\$20.20
No. SF00854	Blue Powder-Free Nitrile Gloves, XL, 100 ea.	\$20.20

Sterile Nitrile Gloves

These long, sterile nitrile gloves are powder-free, minimizing the potential for powder-related complications and come in sterile, sealed pairs.

- The purple color differentiates it from other gloves
- Approved for medical use, 510K
- 9.5 inches long and available small to large



No. SF0084S	Sterile Nitrile Gloves, Sm, 50 sealed pairs	\$98.29
No. SF0084M	Sterile Nitrile Gloves, Md, 50 sealed pairs	\$98.29
No. SF0084L	Sterile Nitrile Gloves, Lg, 50 sealed pairs	\$98.29

Long Cuff Nitrile Gloves

These long cuff, powder-free, durable, 8 mil nitrile gloves provide extra thickness for added protection. Perfect for tough jobs where extra durability or safety is needed. Meets NFPA standards for MES workers and provides protection against bloodborne pathogens. Available in five sizes, and come 50 each per box.



No. SF00871	Long Cuff Nitrile Gloves, Sm, 50 ea.	\$14.44
No. SF00872	Long Cuff Nitrile Gloves, Md, 50 ea.	\$14.44
No. SF00873	Long Cuff Nitrile Gloves, Lg, 50 ea.	\$14.44
No. SF00874	Long Cuff Nitrile Gloves, XL, 50 ea.	\$14.44
No. SF00875	Long Cuff Nitrile Gloves, 2XL, 50 ea.	\$14.44

ResQ-Grip Nitrile Gloves

This is the ultimate in hand protection with superior grip and features our Diamond Pattern Technology. There is no comparison to the super tough nitrile compound that gives ResQ-Grip excellent abrasion resistance and protection. Available in six sizes.



No. SF79720	ResQ-Grip Nitrile Gloves, Sm, 100 ea.	\$20.20
No. SF79730	ResQ-Grip Nitrile Gloves, Md, 100 ea.	\$20.20
No. SF79740	ResQ-Grip Nitrile Gloves, Lg, 100 ea.	\$20.20
No. SF79750	ResQ-Grip Nitrile Gloves, XL, 90 ea.	\$20.20
No. SF79760	ResQ-Grip Nitrile Gloves, 2XL, 90 ea.	\$20.20
No. SF79770	ResQ-Grip Nitrile Gloves, 3XL, 90 ea.	\$20.20

Personal Protection

Latex Gloves

Durable and ambidextrous, latex gloves are a popular choice among public safety professionals and law enforcement.

- Available in four different thicknesses, from 5 mil to 14 mil thick, and two different lengths (10" and 12").
- Available in sizes Small to X-large.



SF0076S	Latex Powder-Free, 10" x 5 mil, Sm, 100 ea.	\$13.81
SF0076M	Latex Powder-Free, 10" x 5 mil, Md, 100 ea.	\$13.81
SF0076	Latex Powder-Free, 10" x 5 mil, Lg, 100 ea.	\$13.81
SF0076XL	Latex Powder-Free, 10" x 5 mil, XL, 100 ea.	\$13.81
SF0082S	Latex Powder-Free, Black, 10" x 6.5 mil, Sm, 100 ea.	\$13.06
SF0082M	Latex Powder-Free, Black, 10" x 6.5 mil, Md, 100 ea.	\$13.06
SF0082L	Latex Powder-Free, Black, 10" x 6.5 mil, Lg, 100 ea.	\$13.06
SF0082XL	Latex Powder-Free, Black, 10" x 6.5 mil, XL, 100 ea.	\$13.06
SF0079S	Latex Powder-Free, Blue, 12" x 10 mil, Sm, 50 ea.	\$20.68
SF0079M	Latex Powder-Free, Blue, 12" x 10 mil, Md, 50 ea.	\$20.68
SF0079	Latex Powder-Free, Blue, 12" x 10 mil, Lg, 50 ea.	\$20.68
SF0079XL	Latex Powder-Free, Blue, 12" x 10 mil, XL, 50 ea.	\$20.68
SF0080S	Latex Powder-Free, 12" x 14 mil, Sm, 50 ea.	\$29.48
SF0080M	Latex Powder-Free, 12" x 14 mil, Md, 50 ea.	\$29.48
SF0080	Latex Powder-Free, 12" x 14 mil, Lg 50 ea.	\$29.48
SF0080XL	Latex Powder-Free, 12" x 14 mil, XL, 50 ea.	\$29.48

SafeGrip Powder-Free Latex Gloves

Four times thicker than standard latex gloves, these gloves provide more protection for high-risk or heavy-duty conditions. Powder free gloves are ideal when handling evidence that you don't want any powder to contaminate. Available in sizes small to XL.



No. SFG375S SafeGrip Powder-Free, Latex Gloves, Sm, 50 ea. \$31.45

No. SFG375M SafeGrip Powder-Free, Latex Gloves, Md, 50 ea. \$31.45

No. SFG375L SafeGrip Powder-Free, Latex Gloves, Lg, 50 ea. \$31.45

No. SFG375XL SafeGrip Powder-Free, Latex Gloves, XL, 50 ea. \$31.45

Heavy-Duty Rubber Gloves

The ideal glove for mixing lab and darkroom chemicals as well as handling strong chemicals and bloody or contaminated evidence at crime scenes. These gloves are 12" long, 18 mil thick, and sold in sets of 12 each.



No. HDG100 Heavy-Duty Rubber Gloves, 12 ea. \$16.11

White Cotton Gloves

In addition to providing protection for fragile evidence, these gloves may be used as liners for latex gloves.

- TIP: Do not wear cotton gloves when using Cyanoacrylate products, serious burns will or can result.
- Available in women's and men's sizes.



No. G4411W White Cotton Gloves, women's, 12 pair \$10.51

No. G4411M White Cotton Gloves, men's, 12 pair \$10.51

Personal Protection

Masks, Respirators, and Facial Protection

Half Mask Particulate Respirators

This is an all-around good mask for investigators and something to always have on hand before entering a crime scene. Different filters can provide different levels of protection:

- Get protection against particulate and airborne pathogens with the High-Efficiency P100 filter.
- Get protection against dust, mists and fumes with the R95 filter.



Often bought with:
P100 Filters, 25 ea.
(SF0034HR),
R95 Filters, 25 ea.
(SF0034DR)



No. SF0034H	Half Mask w/(5 ea.) P100 Filters	\$103.26
No. SF0034D	Half Mask w/(5 ea.) R95 Filters	\$98.29

Multi-Purpose Half Mask Regulator

This versatile mask is comfortable and useful for many different situations with three filter choices. Soft silicone fits around the nose while hard plastic ensures a leak-resistant fit. Filters will protect against almost all decaying matter and hazardous vapors.



No. PU0010	Multi-Purpose Half Mask	\$42.59
No. PU0011	Organic Vapor/Acid Gas Filters, 2 ea.	\$21.27
No. PU0012	P100 Filters (99.7% or airborne particles), 2 ea.	\$17.43
No. PU0013	Organic Vapor/ P100 Filters, 2 ea.	\$37.26

Anti-Putrefaction Mask Kits

Designed to eliminate the overwhelming odor from a decaying body. Investigators can go about the job of investigating, fingerprinting, and photographing – no matter the state of decomposition. Features include:

- Low-maintenance, reusable, and made to fit almost any head size.
- Chemically treated filter, NIOSH/MSHA approved for organic vapors, chlorine, hydrogen chloride, chlorine dioxide, and sulfur dioxide.
- Kits offered with one or two masks



Often bought with:
Organic Vapor/Acid Gas Filters, 2ea.
(PU0011)



No. PU001B	Single Mask Anti-Putrefaction Kit	\$139.38
No. PU003B	Double Mask Anti-Putrefaction Kit	\$265.98

Face Masks

Our multi-layer masks offer high-performance with:

- superior levels of filtration
- comfortable ergonomical shape and nose pads
- climate control valves to help reduce breathing resistance



No. SF0074

No. DRM001

No. DRMV100

No. SF0074	Anti-Odor Mask With Comfort Valve, 1ea.	\$6.13
No. DRM001	Dust Respirator/Mask, 20 ea.	\$31.00
No. DRMV100	Respirator with Comfort Valve, 10 ea.	\$32.73

Personal Protection

Chemical Splash Face Protector

This shield protects the face and eyes from splashed or sprayed chemicals. Shield flips up when not needed.



No. CSFP100 **Chemical Splash Face Protector** **\$40.35**

Chemical Splash Goggles

Protect your eyes from chemicals in the lab. Full coverage goggles fit over regular eye glasses.



No. CSG1001 **Chemical Splash Goggles** **\$4.11**

Jumpsuits and Coveralls

Tyvek® Protection Suits

Each suit is made of Tyvek® fabric, which is resistant to snags, tears, punctures, and many chemicals. Roomy and comfortable.

- Disposable Jumpsuit with Hood. Available in Medium to XX-Large
- Personal Protection Coveralls features shoe covers for head to toe coverage. Available in Large-XX-Large.



Tyvek® Jumpsuits include a hood



Personal Protection Coveralls include a hood and integrated shoe covers

No. SF0072M	Tyvek® Jumpsuit w/hood, Md	\$11.83
No. SF0072	Tyvek® Jumpsuit w/hood, Lg	\$11.83
No. SF0072XL	Tyvek® Jumpsuit w/hood, XL	\$11.83
No. SF0072XXL	Tyvek® Jumpsuit w/hood, 2XL	\$12.90
No. TYV100	Personal Protection Coveralls, Lg	\$10.88
No. TYV101	Personal Protection Coveralls, XL	\$10.88
No. TYV102	Personal Protection Coveralls, 2XL	\$11.94

Tychem® Coverall Jumpsuits

Wear these protective jumpsuits when hazardous materials must be handled.

- Made from DuPont's Tychem® material, it's quality coated with polyethylene, giving it Level B protection against accidental spills, sprays, splashes of acids or bases, and many other chemicals.
- These jumpsuits are also ideal for clandestine drug labs or fire scenes.
- Available in sizes Medium to XX-Large, with hood and attached boots, or without.



No. TYV300M	Coverall w/collar & zipper, M	\$16.47
No. TYV300L	Coverall w/collar & zipper, L	\$16.47
No. TYV300XL	Coverall w/collar & zipper, XL	\$18.60
No. TYV300XXL	Coverall w/collar & zipper, 2XL	\$18.60
No. TYV400M	Coverall w/hood & zipper (boots & elastic wrist), M	\$19.67
No. TYV400L	Coverall w/hood & zipper (boots & elastic wrist), L	\$19.67
No. TYV400XL	Coverall w/hood & zipper (boots & elastic wrist), XL	\$19.67
No. TYV400XXL	Coverall w/hood & zipper (boots & elastic wrist), 2XL	\$21.80

Vented Jumpsuit

Stay cool and comfortable and protect yourself with these jumpsuits. Resistant to snags, tears, and punctures but with breathability and a polypropylene back panel. Available in medium to XX-Large.



No. SF0071M	Vented Jumpsuit, Md	\$17.43
No. SF0071L	Vented Jumpsuit, Lg	\$17.43
No. SF0071XL	Vented Jumpsuit, XL	\$17.43
No. SF0071XXL	Vented Jumpsuit, 2XL	\$17.43

Personal Protection

Micromax® NS Coveralls

Our Micromax® NS Coveralls remain strong, wet or dry and are perfect for work environments with dirt, grime, splashes and spills. They're a general-purpose protective clothing and can be used when hazardous or non-hazardous contaminants are present. Very economical and lightweight.



No. MM100S	Coveralls, Zipper, Sm	\$8.68
No. MM100M	Coveralls, Zipper, Md	\$8.68
No. MM100L	Coveralls, Zipper, L	\$8.68
No. MM100XL	Coveralls, Zipper, XL	\$8.68
No. MM1002XL	Coveralls, Zipper, 2XL	\$8.68
No. MM1003XL	Coveralls, Zipper, 3XL	\$9.86
No. MM1004XL	Coveralls, Zipper, 4XL	\$9.86
No. MM1005XL	Coveralls, Zipper, 5XL	\$9.86
No. MME200S	Coveralls, Zipper, Elastic Wrist & Ankle, Sm	\$9.86
No. MME200M	Coveralls, Zipper, Elastic Wrist & Ankle, Md	\$9.86
No. MME200L	Coveralls, Zipper, Elastic Wrist & Ankle, L	\$9.86
No. MME200XL	Coveralls, Zipper, Elastic Wrist & Ankle, XL	\$9.86
No. MME2002XL	Coveralls, Zipper, Elastic Wrist & Ankle, 2XL	\$9.86
No. MME2003XL	Coveralls, Zipper, Elastic Wrist & Ankle, 3XL	\$10.39
No. MME2004XL	Coveralls, Zipper, Elastic Wrist & Ankle, 4XL	\$10.39
No. MME2005XL	Coveralls, Zipper, Elastic Wrist & Ankle, 5XL	\$10.39
No. MMH300S	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, Sm	\$10.39
No. MMH300M	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, Md	\$10.39
No. MMH300L	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, L	\$10.39
No. MMH300XL	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, XL	\$10.39
No. MMH3002XL	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, 2XL	\$10.39
No. MMH3003XL	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, 3XL	\$10.88
No. MMH3004XL	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, 4XL	\$10.88
No. MMH3005XL	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, 5XL	\$10.88

Shoe Covers

Shoe covers will keep shoes clean and protected and provide extended protection against hazardous materials

Police Embossed Shoe Covers

Protect your shoes with these plastic shoe covers, which feature "POLICE" on the rubber sole. You'll easily know which foot impressions to discount. Available in Medium to X-Large.



No. SF0075M	Embossed Shoe Cover, Md, pair	\$9.81
No. SF0075L	Embossed Shoe Cover, Lg, pair	\$9.81
No. SF0075XL	Embossed Shoe Cover, XL, pair	\$10.88

Disposable Tyvek® Shoe Covers

One size fits all and keeps shoes clean and protected.



No. SF0073	Disposable Shoe Covers, pair	\$2.72
------------	------------------------------	--------

Disposable Protective High-Top Boots

These Tyvek® boots have vinyl, non-skid soles. Available in medium or large.



No. TYV200	High-Top Boots, Md, pair	\$4.85
No. TYV201	High-Top Boots, Lg, pair	\$4.85

Micromax® NS Shoe & Boot Covers

Breathable for comfort, they remain strong whether they're wet or dry.



No. SC100SM	Shoe Covers, Elastic Ankle, Sm/Md, pair	\$1.66
No. SC100LXL	Shoe Covers, Elastic Ankle, L/XL, pair	\$1.66
No. SC1002XL	Shoe Covers, Elastic Ankle, 2XL, pair	\$1.66
No. BC17SM	Boot Covers, Elastic Top, Sm/Md, pair	\$1.95
No. BC17LXL	Boot Covers, Elastic Top, L/XL, pair	\$1.95
No. BC172XL	Boot Covers, Elastic Top, 2XL, pair	\$2.02

NARKsafe™ Protection & Packaging

NARK Single-Use Powder Substance Packaging, Transport, Testing Kit

Any powdery substance encountered on the street should be considered dangerous. Law enforcement needs to proceed with caution while handling, transporting, and testing these substances.

Following DEA and CDC guidelines, we developed this one-time use PPE Kit to be used specifically around powdery substances that could contain fentanyl, carfentanil, or other powerful opioid-based substances.

Kit contains: coveralls; P-100 respirator; two pair of nitrile gloves; goggles; evidence bag; a biohazard bag; and donning/doffing instruction card.

No. NARK300PPE **NARK Single-Use Powder Substance Kit** **\$31.00**



NARKsafe Patrol Kit

NARKsafe™ Patrol Kit keeps you safe when you encounter narcotics in the field. The kit includes one respirator mask with 100% filtration of harmful substances, fully sealed goggles and a pair of nitrile gloves. To stay safe, take the same precautions on patrol as they do in the lab, whether conducting presumptive field tests or collecting drug evidence, wear protective gear. *If fentanyl is suspected, use our NARK300PPE kit for more complete protection.*

No. NARK300PK **NARKsafe Patrol Kit** **\$22.95**



Disposable P100 Respirator Mask

This Respirator offers the highest level of NIOSH rated filtration efficiency in a maintenance-free respirator. It features a 3M™ Cool Flow™ Exhalation Valve, a facesal, adjustable noseclip, lightweight construction and an all filter design. Sold individually. Meets 42 CFR84 standard.

No. NARK300P100R **Disposable P100 Respirator Mask** **\$10.30**



Chemical Splash Goggles

Protect your eyes from chemicals in the lab. Full coverage goggles fit over regular eye glasses.



No. CSG1001 **Chemical Splash Goggles** **\$4.11**

Black Powder-Free Nitrile Gloves

These comfortable, ambidextrous nitrile gloves provide extra puncture and abrasion resistance. Always use nitrile gloves while working with illicit substances and NARK presumptive drug kits.



No. SF0081S	Black, Powder-Free, Sm, 100 ea.	\$20.20
No. SF0081M	Black, Powder-Free, Md, 100 ea.	\$20.20
No. SF0081L	Black, Powder-Free, Lg, 100 ea.	\$20.20
No. SF0081XL	Black, Powder-Free, XL, 100 ea.	\$20.20

Long Cuff Nitrile Gloves

These long cuff, powder-free, durable, 8 mil nitrile gloves provide extra thickness for added protection. Perfect for tough jobs where extra durability or safety is needed. Meets NFPA standards for MES workers and provides protection against bloodborne pathogens.



No. SF00871	Long Cuff Nitrile Gloves, Sm, 50 ea.	\$14.44
No. SF00872	Long Cuff Nitrile Gloves, Md, 50 ea.	\$14.44
No. SF00873	Long Cuff Nitrile Gloves, Lg, 50 ea.	\$14.44
No. SF00874	Long Cuff Nitrile Gloves, XL, 50 ea.	\$14.44
No. SF00875	Long Cuff Nitrile Gloves, 2XL, 50 ea.	\$14.44

Micromax® NS Boot Covers

Breathable for comfort, they remain strong whether they're wet or dry.



No. BC17SM	Boot Covers, Elastic Top, Sm/Md, pair	\$2.02
No. BC17LXL	Boot Covers, Elastic Top, L/XL, pair	\$2.02
No. BC172XL	Boot Covers, Elastic Top, 2XL, pair	\$2.02

NARKsafe™ Protection, Packaging & Disposal

Micromax® NS Coveralls

According to the Drug Enforcement Agency (DEA) and Centers for Disease Control (CDC), a full coverall with hood is recommended for handling dangerous illicit substances. Below are our most popular sizes. You can find our full size availability on page 8.




Often bought with:
Micromax NS
Shoe & Boot Covers



No. MMH300M	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, Md	\$10.04
No. MMH300L	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, L	\$10.04
No. MMH300XL	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, XL	\$10.04
No. MMH3002XL	Coveralls, Zipper, Attached Hood, Elastic Wrist & Ankle, 2XL	\$10.04

Disposable Hazardous Material Bags

Properly dispose of your used personal PPE in our 18" x 30 1/2" hazardous material bags: easily identifiable by its bright yellow color. Always dispose of any hazardous waste according to Federal, state, and local regulations. Sold in packs of 10 bags each.

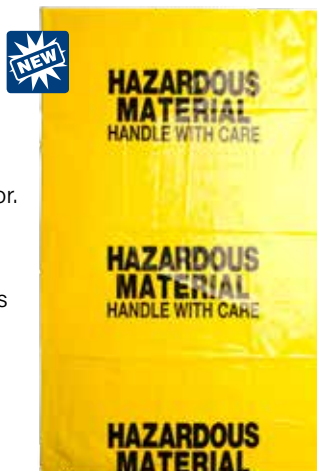

No. NARK300HVB	Disposable Hazard Bags, 18" x 30 1/2", 10 ea.	\$22.95
----------------	-----------------------------------------------	---------

Fentanyl Warning Labels

These attention-getting labels are printed in full color on durable crack-n-peel material, and measure 2.875" x 3.75".



No. NARK300FENL	Fentanyl Warning Labels, 100/pk	\$5.95
-----------------	---------------------------------	--------

Keepsafe Poly Evidence Bags

These bags feature a high degree of security: once closed, unauthorized attempts to break the seal will cause "STOP" signs to appear in the orange strips of the tape. Also includes a tear-away tab for evidence tracking.



No. KSB610C	Keepsafe Bags, 6" x 10 1/4", 100 ea.	\$42.59
No. KSB1013C	Keepsafe Bags, 10" x 13", 100 ea.	\$53.25
No. KSB1216C	Keepsafe Bags, 12 x 16 7/8", 100 ea.	\$63.91

Write-Block Zip-Top Evidence Bag

These zip-top bags are made of 0.004" thick polyethylene and feature a write-on block for easy labeling of contents. Sold in packs of 100 each.



No. ZWB02	Zip-Top Write Block Evidence Bag, 9" x 12", 100 ea.	\$27.67
-----------	-----------------------------------------------------	---------

Sharpstube™ Syringe Transport Tubes

Designed to hold a hypodermic syringe for transportation, these tubes are well-marked with bright orange-red international biohazard markings. Sold in pack of 12 each.



No. ECT2	Syringe Transport Tube, 12 ea.	\$23.75
----------	--------------------------------	---------

NARKsafe NDA Narcotic Disposal Agent

NARKsafe NDA is a drug disposal solution that enables the onsite disposal of drug evidence (20 grams or less) – quickly, easily and safely, rendering the drugs unusable.

The only polymer-based drug disposal solution on the market, NARKsafe NDA enables onsite disposal in an effective, environmentally-safe method that can be tossed into the trash.



No. NDAK100	NARKsafe NDA Kit, 1 each	\$10.85
-------------	--------------------------	---------

NARKsafe™ Protection & Packaging

Liquid Sampling Kits

Many protocols for the collection, transport, and storage of hazardous substances, including those often found at clandestine meth labs, recommend or require placing the substance in a glass jar and then placing the jar inside a Nalgene® jar. The use of leak proof glass and Nalgene® jars helps to enhance the safety of all personnel and preserve the evidence. Each of the Sirchie Liquid Sampling Kits is designed to have a glass jar that fits perfectly into a Nalgene® jar with pipettes for collection of the substance.

30ml Liquid Sampling Kit

This kit includes 30ml Nalgene® Bottles (12 each), 4ml Glass Vials (12 each), and Disposable Pipettes (12 each).



No. LSK100 **30ml Liquid Sampling Kit** **\$45.50**

60ml Liquid Sampling Kit

This kit includes 60ml Nalgene® Bottles (12 each), 12ml Glass Vials (12 each), and Disposable Pipettes (12 each).



No. LSK200 **60ml Liquid Sampling Kit** **\$53.95**

125ml Liquid Sampling Kit

This kit includes 125ml Nalgene® Bottles (12 each), 16ml Glass Vials (12 each), and Disposable Pipettes (12 each).



No. LSK300 **125ml Liquid Sampling Kit** **\$54.95**

250ml Liquid Sampling Kit

This kit includes 250ml Nalgene® Bottles (12 each), 30ml Glass Vials (12 each), and Disposable Pipettes (12 each).



No. LSK400 **250ml Liquid Sampling Kit** **\$69.95**

Disposable Liquid Sampling Pipettes

Use these disposable pipettes for sampling unknown liquids. Store collected liquids in glass vials placed inside Nalgene bottles.



No. PIP25 **Disposable Liquid Sampling Pipettes, 25 ea.** **\$2.50**

Translucent Nalgene® Bottles

You have the option to get the wide-mouth Nalgene® bottles separately.



No. NAJAR0001 **Translucent Nalgene® Bottles, 30ml, 12 ea.** **\$16.39**

No. NAJAR0002 **Translucent Nalgene® Bottles, 60ml, 12 ea.** **\$18.04**

No. NAJAR0004 **Translucent Nalgene® Bottles, 125ml, 12 ea.** **\$20.63**

No. NAJAR0008 **Translucent Nalgene® Bottles, 250ml, 12 ea.** **\$28.89**

No. NAJAR0016 **Translucent Nalgene® Bottles, 500ml, 12 ea.** **\$42.61**

No. NAJAR0032 **Translucent Nalgene® Bottles, 1 liter, 6 ea.** **\$42.55**

Amber Nalgene® Lab Quality Bottles

Store, ship and package light-sensitive liquid using these narrow-mouth Amber HDPE lab quality Nalgene® bottles. Featuring excellent chemical resistance to most acids, bases, and alcohols.



No. NABOT9125 **Amber Nalgene® Bottles, 4ml, 332 ea.** **\$283.52**

No. NABOT9025 **Amber Nalgene® Bottles, 8ml, 12 ea.** **\$21.59**

No. NABOT9050 **Amber Nalgene® Bottles, 15ml, 150 ea.** **\$133.53**

No. NABOT0001 **Amber Nalgene® Bottles, 30ml, 12 ea.** **\$42.91**

No. NABOT0002 **Amber Nalgene® Bottles, 60ml, 12 ea.** **\$45.31**

No. NABOT0004 **Amber Nalgene® Bottles, 125ml, 12 ea.** **\$54.64**

No. NABOT0008 **Amber Nalgene® Bottles, 250ml, 12 ea.** **\$73.83**

No. NABOT0016 **Amber Nalgene® Bottles, 500ml, 12 ea.** **\$109.01**

No. NABOT0032 **Amber Nalgene® Bottles, 1 liter, 6 ea.** **\$91.63**


No. NABOT0010 **Amber Nalgene® Bottles, 1 gal, 1 ea.** **\$48.99**

NARKsafe™ Testing Workstation



NARKsafe™ Illicit Substance Testing Workstation

Fentanyl-related deaths are increasing as the opioid epidemic increases worldwide. Fentanyl is deadly if swallowed, inhaled, or absorbed through the skin, so being protected when handling suspected substances is paramount. NARKsafe™ fume hood enclosures are ideal for protecting professionals during the forensic analysis of fentanyl and other dangerous substances. Only this system offers a turbulence-free, horizontal laminar airflow with dual HEPA filtration.

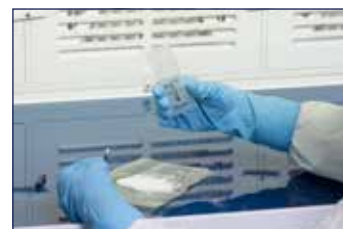


NARKsafe™ enclosure provides a negative pressure environment to pull any particulate away from the operator and eliminate end-user exposure. It also allows for the addition of a bonded carbon filter, so analysis of other chemical solutions or marijuana exhibits can be safely handled in the same workspace. Available in three widths: 36", 48", 72".

Why NARKsafe™ is the only enclosure to use when testing illicit substances:

- Horizontal laminar airflow pattern
- Dual HEPA filtration
- HEPA filters with 99.997% efficiency rating for collecting 0.3-micron particles
- Electrostatically charged pre-filter to collect larger particulate and minimize clogging of the HEPA filter
- ASHRAE 110 testing adjustable face velocity
- Ability to replace contaminated filters under negative pressure.
- Blue flooring to see substances clearly when testing



i When testing illicit substances and changing filters, the operator should always follow Drug Enforcement Administration (DEA) and the Centers for Disease Control (CDC) guidelines for handling fentanyl, carfentanil and acetyl-fentanyl.

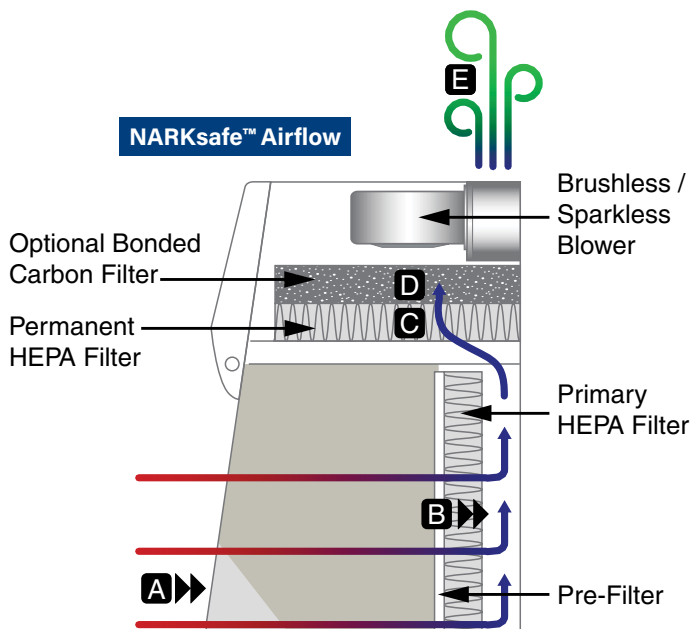

No. ACDD760	NARKsafe™ 36" Powder Safety Workstation, 110V AC w/prefilter, carbon filter & HEPA filters	\$9,740.00
No. ACDD760220	NARKsafe™ 36" Powder Safety Workstation, 220V AC w/prefilter, carbon filter & HEPA filters	\$9,740.00
No. ACDD770	NARKsafe™ 48" Powder Safety Workstation, 110V AC w/prefilter, carbon filter & HEPA filters	\$10,650.00
No. ACDD770220	NARKsafe™ 48" Powder Safety Workstation, 220V AC w/prefilter, carbon filter & HEPA filters	\$10,650.00
No. ACDD780	NARKsafe™ 72" Powder Safety Workstation, 110V AC w/prefilter, carbon filter & HEPA filters	\$14,575.00
No. ACDD780220	NARKsafe™ 72" Powder Safety Workstation, 220V AC w/prefilter, carbon filter & HEPA filters	\$14,575.00

Presumptive Drug Identification

NARK II Presumptive Drug Testing Reagent Pouches

Our NARK drug identification pouches can presumptively identify several families of substances that could be abused drugs. Each of the tests is comprised of one or more chemical reagents. When a color or series of colors occurs within a specific testing sequence, a positive identification may be presumed.



Most popular reagent tests:

No.	Description	Presumptively tests for	Price
NARK2005	Duquenois-Levine Reagent, 10/box	Marijuana, Hash, Hash Oil, THC Concentrates	\$21.95
NARK2007	Scott Reagent (Modified), 10/box	Cocaine Salts (HCl), Cocaine Base (Crack)	\$21.95
NARK20011	Mecke's Modified Reagent, 10/box	Heroin (White, Brown, Black Tar)	\$21.95
NARK20015	Methamphetamine/MDMA Reagent, 10/box	Methamphetamine, MDMA (Ecstasy)	\$21.95
NARK20023	Synthetic Cannabinoids Reagent, 5/box	Synthetic Cannabinoids (indole formulations)	\$21.95
NARK20033	Fentanyl Reagent, 10/box	Fentanyl, Acetyl-fentanyl, Heroin	\$21.95

No.	Description	Presumptively tests for	Price
NARK2001	Marquis Reagent, 10/box	Opium Alkaloids (Heroin/Morphine), MDMA (Ecstasy), Amphetamines/Meth	\$17.00
NARK2002	Nitric Acid Reagent, 10/box	Differentiate Heroin/Morphine	\$17.00
NARK2003	Dille-Koppanyi Reagent, 10box	Barbiturates	\$21.95
NARK2004	Ehrlich's Reagent, 10/box	LSD	\$21.95
NARK2008	Methadone Reagent, 10/box	Methadone	\$21.95
NARK2009	PCP/Methaqualone Reagent, 10/box	PCP, Methaqualone	\$21.95
NARK20010	Special Opiates Reagent, 10/box	Codeine, Heroin, Morphine, Oxycodone, Buprenorphine	\$21.95
NARK20012	Froehde's Reagent, 10/box	Talwin® (Pentazocine), DMT (Dimethyltryptamine)	\$21.95
NARK20013	Ephedrine Reagent, 10/box	Ephedrine, Pseudoephedrine	\$21.95
NARK20014	Valium Reagent, 10/box	Valium, Rohypnol, Ketamine	\$21.95
NARK20019	Mayer's Reagent, 10/box	General Screening, Narcotic Alkaloids	\$17.00
NARK20020	KN Reagent (Fast Blue B Salts), 10/box	Marijuana, Hashish, Hash Oil, THC Concentrates (DAB, WAX, BHO)	\$21.95
NARK20021	GHB Reagent, 5/box	GHB (Gamma Hydroxybutyrate)	\$21.95
NARK20022	Mandelin Reagent, 10/box	General screening for Methadone, Amphetamine-Type Compounds	\$17.00
NARK20024	MDPV Reagent, 10/box	MDPV, Butylone, Ethylone (bk-MDEA), Methylone, Naphyrone, Pentylone	\$21.95
NARK20025	Mephedrone Reagent, 10/box	Mephedrone (Synthetic Cathinones)	\$21.95
NARK20026	a-PVP Reagent, 10/box	a-PVP (Flakka)	\$21.95
NARK20029	2C Reagent, 10/box	2C-B, 2C-B-NBOMe, 2C-C, 2C-C-NBOMe, 2C-E, 2C-E-NBOMe, 2C-I, 2C-I-NBOMe, 2C-N, 2C-T-7	\$21.95
NARK20030	Psilocybin/Psilocin Reagent, 10/box	Psilocybin, Psilocin (Mushrooms)	\$21.95
NARK20031	Liebermann Reagent, 10/box	Methamphetamine, Morphine	\$17.00
NARK20032	Mollies Reagent, 10/box	Mollies	\$17.00
NARK20033	Fentanyl Reagent, 10/box	Fentanyl, Acetyl-fentanyl, Heroin	\$21.95

No.	Description	Price
NARK2006	Acid Neutralizer Liquid, 1 oz bottle	\$5.95
NARK200N	Acid Neutralizer Powder, 0.5 oz. pouch, 10 ea.	\$11.95

Presumptive Drug Identification

NARK Test Reagents and Components

These reagents are in disposable tube format with each reagent in a glass ampoule.

Our top 3 drug reagents:

- NARK Meckes Modified Reagent (Heroin.) No. NAR10016
- NARK Special Opiates (Codeine, Heroin, Oxycodone). No. NAR10022
- NARK MDPV Reagent (Synthetic Cathinones). No. NAR10024



No.	Description	Presumptively tests for	Price
NAR10001	Mayer's Reagent, 10/box	Narcotic Alkaloids	\$15.95
NAR10002	Marquis Reagent, 10/box	Opium Alkaloids	\$15.95
NAR10003	Nitric Acid Reagent, 10/box	Differentiate Heroin/Morphine	\$15.95
NAR10004	Cobalt-Thiocyanate Reagent, 10/box	Cocaine/Crack	\$15.95
NAR10005	Dille-Koppanyi Reagent, 10/box	Barbiturates	\$15.95
NAR10006	Mandelin Reagent, 10/box	Amphetamines	\$15.95
NAR10007	Ehrlich's Reagent, 10/box	LSD	\$15.95
NAR10008	Duquenois Reagent, 10/box	Marijuana, Hash, Hash Oil, THC Concentrates	\$15.95
NAR10009	KN Reagent, 10/box	Marijuana (green plant & seeds)	\$15.95
NAR10013	Modified Cobalt-Thiocyanate Reagent, 10/box	Cocaine Freebase, Crack	\$15.95
NAR10014	Valium Reagent, 10/box	Valium, Rohypnol, Methcathinone, Ketamine	\$15.95
NAR10015	Meth/MDMA Reagent, 10/box	Methamphetamine, MDMA (Ecstasy)	\$15.95
 NAR10016	Mecke's Modified Reagent, 10/box	Heroin (white, brown, black tar)	\$15.95
NAR10018	Froehde's Reagent, 10/box	Talwin® (Pentazocine), DMT (Dimethyltryptamine)	\$15.95
NAR10019	Ephedrine Reagent, 10/box	Ephedrine, Pseudoephedrine	\$15.95
 NAR10022	Special Opiates Reagent, 10/box	Codeine, Heroin, Oxycodone	\$15.95
 NAR10024	MDPV Reagent, 10/box	MDPV (Synthetic Cathinones)	\$15.95
NAR10025	Mephedrone Reagent, 10/box	Mephedrone (Synthetic Cathinones)	\$15.95
NAR10028	GHB Reagent, 5/box	GHB (Gamma Hydroxybutyrate)	\$15.95
NAR10029	PCP/Methaqualone Reagent, 10/box	PCP, Methaqualone	\$15.95
NAR10030	Psilocybin/Psilocin Reagent, 10/box	Psilocybin, Psilocin (mushrooms)	\$15.95
NAR10010	Polyester Bags, 10/box		\$2.95
NAR10011	Acid Neutralizer, 1 oz. Bottle		\$5.95

NARK Analysis Reagent Kit

We included reagents for the most frequently abused drugs in this field test kit so you always have them on hand. Each individual reagent is sealed in a thin, glass ampoule which eliminates cross-contamination with other reagents.

Reagents included for: Marijuana, GHB, Hashish, THC, Cocaine, Crack/Freebase, Amphetamines, Barbiturates, Narcotic Alkaloids, Opium Alkaloids, PCP, Hallucinogens, Heroin, Morphine.

No. NAR100 **NARK Analysis Reagent Kit** **\$134.95**



Alternate Light Sources

megaMAXX ALS System

This system has everything an investigator needs at the scene or in the lab for finding evidence with ALS. The kit holds seven self-contained light sources. Other components are:

- Light Diffuser. A softer light is required for oblique viewing and photography. This diffuser will soften the LED lights
- ForensiTorch UV Light. This light should always be used in preliminary examining and locating evidence
- megaMAXX Tripod. You need both hands free during an investigation to take photos or process an area. This tripod has a flexible extension and can be precisely positioned.



No. MMX300

megaMAXX ALS System

\$1,439.17

Tactical Max Forensic Light Kit

The Tactical Max provides all the light needed to conduct an on-scene forensic investigation: white, 365nm UV, and 455nm blue. Each handheld forensic light is made of tough, anodized aluminum and contains an extremely bright 1-Watt LED. This kit was built for field use, contained in a heavy duty nylon case and includes a flexible tripod that can be used as a third hand. Lightweight and portable, yet rugged and durable—there's no reason to not take the Tactical Max Forensic Light Kit everywhere. The Tactical Max Flexible Mini Tripod can be purchased separately and is compatible with our popular MMX300 series forensic lights!



No. TMX200

Tactical Max Forensic Light Kit

\$388.07

No. TMXFT1

Flexible Tripod (only)

\$54.59

Forensic Alternate Light Source

Physical evidence is the backbone of criminal investigations. If the perpetrator has left some form of evidence behind, this tool can help you find it, tie a potential perpetrator to the scene, and help prosecutors file formal charges.



Features include:

- Fully portable for unlimited movement around a crime scene. It can be hand-held or mounted on a tripod and operates by battery.
- Filters included with the light source have a selection of six different wavelengths, covering 365nm to 530nm.
- A high-intensity xenon lamp which is ideal to search for hair,



fibers, footwear evidence, and more.

No. FAL3000	Forensic Alternate Light Source	\$4,258.87
No. FAL203	Replacement High-Intensity Xenon Lamp	\$224.01
No. FAL306	AC Adapter Kit 120V/60Hz	\$273.18
No. FAL306220	AC Adapter Kit 220V/50Hz	\$273.18
No. FAL309	Rechargeable 12V DC BAh Battery Pack	\$228.92
No. FAL309B	Replacement 12V DC BAh Battery Set	\$131.07

UV ForensiTorch

If investigators conduct a search in subdued light or broad daylight, the ForensiTorch is ideal for the preliminary examination and location of forensic evidence. It's a single, very high intensity 365nm longwave UV LED and has 50,000 hours of life. **Improved design includes a USB cable for charging the included lithium ion rechargeable battery.**



No. UVFT100A

UV ForensiTorch

\$92.12

Evidence Photography

Investigators and witnesses can describe scenes, sometimes in great detail. However, photographs can tell the story better with objectivity. They freeze time and leave a permanent record of the scene. We've put together several options of photography kits to meet the needs of your department.

Evidence Photography Kits

We offer three levels of photography kits to meet the needs of different agencies and departments. They include:

Basic Photography Kit

A point-and-shoot camera with a variety of forensic scales. Photographing crime scenes is so important, all departments should have this standard kit.



No. CAM100

Basic Photography Kit

[Call for Price](#)

Advanced Photography Kit

An advanced compact digital camera with a variety of forensic scales for proper documentation. This camera gives more seasoned investigators more functions, such as an LCD screen and an ability to shoot video.



No. CAM200

Advanced Photography Kit

[Call for Price](#)

Professional Photography Kit

We've put together the latest in technology to provide CSI with all the photographic equipment they'll need to properly and efficiently photograph a scene. This kit also covers all visible light applications, from overall scene captures to life-sized macrophotography. With the TTL flash cord (included), you'll be able to use the speedlight flash (also included) off-camera for oblique lighting.



No. CAM300

Professional Photography Kit

[Call for Price](#)

Comprehensive Crime Scene Photography



You will leave this class knowing:

- How to take proper crime scene photographs in a variety of environments
- How to present crime scene photos in judicial and administrative proceedings
- How to train other professionals in taking photographs



For more information about this and all other training courses, visit www.sirchie.com or call us at (800) 356-7311.





Evidence Photography

SirchStand Forensic Photography Stand

One of our more popular products, we designed the SirchStand to make forensic photographers' jobs easier. It's an easy-to-use stand with a camera mount suitable for most cameras. It includes flexible arms installed on the base to hold additional forensic lights.

The SirchStand features:

- Easy to adjust column for one-hand operation to raise or lower the camera.
- Sturdy column that can support most mid-range to large cameras
- A magnetic base with a 1 cm square grid for scaling appropriately
- Flexible arms on the base and the camera mount. Our TMX LED forensic lights (TMX200) can be used for oblique light photography. UV 365nm and 455nm blue lights can be used for photographing bodily fluid stains, enhanced cyanoacrylate fuming results, or fluorescent powder developed latent prints
- Portable, weighing only 22 lbs with a removable column.

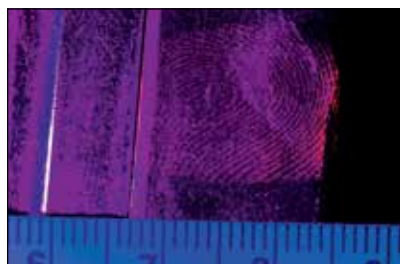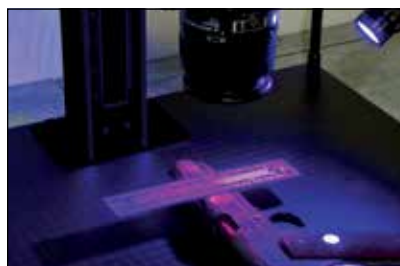

No. PHTSS1000A

SirchStand Forensic Photography Stand

\$1,547.33

Latent Fingerprint Development

Latent prints, left at a crime scene, may or may not be visible to the naked eye. Prints can provide recognizable ridge details that can help identify the owner of those prints. The print residue that is left behind can often be very fragile so caution should be taken to select and execute latent print processing methods that don't damage any possible ridge detail. As with any evidence, all visible print detail should be photographed immediately.

Overview of Chemical Tools and Methods

It's essential that technicians understand the sequence of development products since developmental reactions will change the chemical nature of the fingerprint and it might become impossible to develop.

Method	Why it's useful to you
Powders	Powders adhere to water and fatty deposits and generally work better with newer fingerprints. Choose a powder that contrasts with the background.
Fluorescent Powders	Fluorescent powders are especially useful on patterned backgrounds. An ultraviolet or forensic light source and filter goggles are recommended.
Small Particle Reagent	It adheres to the fatty components of skin secretions. The reaction product is very fragile and should be photographed immediately.
Cyanoacrylate (glue fuming)	This method produces excellent results on Styrofoam and plastic bags. It reacts with moisture to form a hard, whitish deposit.
Adhesive Side Powder	Adhesive side powder attaches itself to print impressions on the adhesive side of tapes. It produces excellent results on light-colored or transparent tapes including duct tapes, masking tapes, plastic surgical tape, foam tape, clear and reinforced packing tapes and labels.

Guide to Surfaces

Surface	Treatment Method
Smooth, non-porous: glass, hard plastic mouldings, surfaces treated with paint or varnish	Powders, iodine, small particle reagent, cyanoacrylate/fluorescent dyes
Rough, non-porous: rough, textured surfaces, grained plastic mouldings	Small particle reagent, cyanoacrylate/fluorescent dyes
Paper & cardboard: paper & cardboard that hasn't been waxed or plastic-coated.	Iodine, ninhydrin, DFO, silver nitrate, or physical developer. Powders won't work with older fingerprints
Plastic packaging media: Polyethylene, polypropylene, cellulose acetate, and laminated paper surfaces	Iodine, small particle reagent, cyanoacrylate/fluorescent dyes, and powders
Soft vinyl (PVC), rubber and leather: simulated leather and cling film	Iodine, small particle reagent, cyanoacrylate, and powders
Untreated metal: untreated, bare metal surfaces; metal surfaces that haven't been painted or lacquered.	Small particle reagent, powders, cyanoacrylate / fluorescent dyes and powders
Unfinished wood: wood surfaces that haven't been painted or treated.	Ninhydrin. You can use powders on smooth wood and silver nitrate or physical developer on light woods.
Wax and waxed surfaces: items made of wax, such as candles, wax-coated paper, cardboard, wood surfaces	Nonmetallic powders, cyanoacrylate / fluorescent dyes
Adhesive-coated surfaces: tapes and similar surfaces that are not likely to dissolve in water	Adhesive-side powders

Latent Fingerprint Powders and Applicators

Our Fingerprint Series

Moving through our Fingerprint classes will give you the knowledge and confidence to classify fingerprints, compare fingerprints and palm prints, and ultimately, pass the IAI Latent Print Exam.



- Introduction to the Science of Fingerprints
- Basic Latent Fingerprint Comparison
- Comprehensive Advanced Latent Print Comparison
- Advanced Latent Ridgeology
- Advanced Latent Palm Print Comparison
- Mastering the IAI Latent Print Exam



Regular Latent Print Powders

Our print powders offer the best sensitivity, adherence qualities, and reproduction capabilities of any powders available. Use them to produce excellent results with many color options to choose from.



Latent Powder Brushes

We have an extensive line of brushes to meet your needs, whether it's tight spaces or dusting large areas. Our brush line includes Carbosmoove, Goat Hair, Max-Spread, Fiberglass, Fiberglass Mini Light, and Feather Dusters. Contact your Sirchie Sales Representative and they'll be able to help you find the right brush for the job.



Magnetic Latent Print Powders

Our high-quality magnetic powders revolutionized the way crime scene technicians approached the task of developing latent prints. Easy-to-use, magnetic powder has become a popular choice among top investigators.



Magnetic Wands

One main advantage of using magnetic wands with magnetic powder is that it's easy to clear surfaces of excess powder, giving maximum photographic effectiveness. Our magnetic applicators can be used with any magnetic powder.



Latent Fingerprint Powders

DNA-Free Latent Print Processing

With advances in DNA analysis, full profiles can be obtained from very small samples. This has led to CSI processing crime scenes for "touch DNA". This means DNA that is recovered from skin cells that are left behind when a person touches or comes into contact with something, such as clothes or weapons.

It is therefore critical that the integrity of the crime scene be maintained, including eliminating cross-contamination possibilities from the tools and methods used to collect the DNA and the tools used to process the crime scenes.



To meet your needs, we now offer DNA-Free Processing tools, including brushes, powders, and lifting tape. They operate and have the same integrity as our products but:

- These items are treated with a scientifically proven DNA destroying process that penetrates and decontaminates throughout, not just on the surface.
- Each lot is certified through a third party testing to ensure each batch processed is DNA free.
- All items are packaged to prevent contamination before use, and are meant to be used only once.

DNA-Free Silk Black Latent Print Powder

This powder offers superior contrast on light-colored surfaces.

No. DNA101L Silk Black, 1 oz. \$6.95

DNA-Free Indestructible White Latent Print Powder

Indestructible White is the whitest of white powders available from any source.

No. DNA103L Indestructible White, 1 oz. \$6.95

DNA-Free Silver/Black Latent Print Powder

Silver/Black latent print powder is a two-color powder that will contrast with any color surface.

No. DNASB201L Silver/Black, 1 oz. \$6.95

DNA-Free Fiberglass Brush

These fiberglass brushes are the same as our popular 122L brush, but packaged in a breathable pouch, and processed and certified to be DNA free.

No. DNA122L Fiberglass Brush \$7.41

DNA-Free Squirrel Hair Brush

These squirrel hair brushes are the same as our 120LS brush, but packaged in a breathable pouch, and processed and certified to be DNA free.

No. DNA120LS Squirrel Hair Brush \$21.31


DNA-Free Lifting Tape

Pressure wound, heat treated and "fish-eye" free, this lifting tape is also treated to destroy DNA and is certified DNA free.

No. DNA144L2 Lifting Tape, 2" x 120", 1 roll \$9.54

DNA-Free Latent Print Kits

This single use kit has everything you need to process, lift, and preserve fingerprints at the crime scene while ensuring that no outside DNA is introduced. Kits are processed and certified as DNA Free.

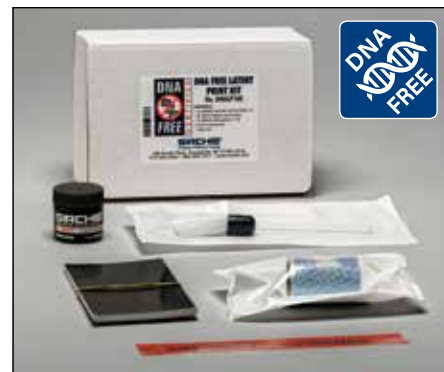

No. DNALP100 w/Fiberglass Brush \$21.80

No. DNALP100S w/Squirrel Hair Brush \$32.73



For more information on DNA-Free products call us at 800.356.7311 or visit www.sirchie.com

RUVIS Systems


Reflected Ultraviolet Imaging

Detect latent fingerprints without the use of powders or chemicals by using Reflected Ultraviolet Imaging Systems (RUVIS) at the scene or in the lab. RUVIS is most effective on smooth, non-porous surfaces (flooring, walls, countertops, tables) and on multi-colored surfaces like magazine covers.

How it works: RUVIS uses shortwave UV light. This light is reflected or absorbed by most surfaces but with the residues in fingerprints, the UV light is reflected and scattered. Using a special lens and a UV filter, the system converts the UV wavelength to visible light, revealing the latent print.

You need RUVIS to:

- Find latent fingerprints and footwear impressions on smooth, non-porous surfaces
- Enhance the ability to see cyanoacrylate fumed prints without further enhancements.
- Detect other evidence you may not have seen otherwise.

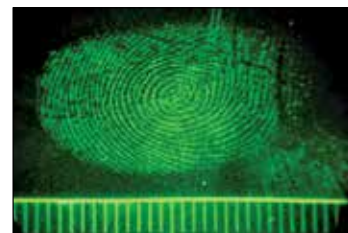

RUVIS Pro Tips



Smudges and smears are commonplace at the crime scene and are often passed over by evidence technicians due to lack of identifiable ridge detail. However, deposits and residues can still be useful if they contain the perpetrator's DNA profile. There's been a number of cases when this has helped to identify evidence.

For evidence technicians, the problem is that you have to locate these potentially flawed latent prints. The KRIMESITE Imager can find these untreated latent fingerprints on many smooth non-porous surfaces.

Please keep in mind that DNA will degrade if DNA is overexposed to shortwave UV light for an extended period of time. The degradation depends on many factors, including:

- The intensity of the light source being used
- The distance that the light source is from the DNA
- The amount of exposure time
- The type of stain or residue.

Your agency should develop a standard operating procedure (SOP), such as the one below:

If you're using the Imager to search for latent prints and see what may be biological stains or other sources of trace DNA, we highly recommend you follow the procedure below:

- Turn OFF Shortwave UV light source
- Swab the suspected area for potential trace DNA
- Cover suspected area with UV-blocking material to shield potential trace DNA from shortwave UV light source.
- Finish conducting your search for untreated latent fingerprints.



Training is essential to get the most out of your RUVIS unit. Sirchie offers training for purchase as a service via our KSSTRAIN offering. Contact us for more details.

RUVIS Systems

KRIMESITE™ Direct View Kit

From its lightweight body to its superior lens with dual filter slide assembly, the KRIMESITE™ Direct View sets a new standard for crime scene search capability.

The 60mm UV lens is our very finest and doesn't require any add-on lens adapters for close-up or distance viewing. Because of this exceptional lens and filter combination, you

can view close-up and distance and photograph without ever changing the mechanical configuration of the device. This kit comes with everything you need to search a crime scene for latent prints.



The KRIMESITE™ Imager can be configured for use several ways.

No. KSS60 **KRIMESITE™ Direct View Kit** **\$15,007.50**

KRIMESITE™ Direct View Kit with Black Talon Stabilizer

The Black Talon Stabilizer is designed for investigators to get the correct lens-to-subject distance for photography. Available as a single item or as part of the KRIMESITE™ Imager Direct View Kit.



Features are:

- It has a built-in precision laser-engraved sliding photo evidence scale that is easily positioned into or out of the field-of-view, allowing the photography to be scaled 1-to-1 at a later time.
- Dual 4-watt shortwave UV lamps are attached to the fully adjustable support/swing arms so you can get the angle you need for maximum contrast.
- The talon stabilizer was created for the KRIMESITE™ Imager and it can be used as a rest spot for the Imager when not in use.
- For close-up photography, tripods can become cumbersome. The stabilizer eliminates this with sturdy support on vertical or horizontal surfaces/evidence.
- For longer shots when a standard tripod is needed, the stabilizer can be quickly mounted to the tripod.



No. KSS60TALON **Direct View Kit w/Black Talon Stabilizer** **\$15,473.25**
No. BTS100 **Black Talon Stabilizer (only)** **\$491.72**

RUVIS Systems

Scan-n-Find™ Pocket Imager Kit

An affordable alternative to the KRIMESITE™ Imager, the Scan-n-Find Imager simplifies the search for potential latent print evidence that could be easily overlooked.

- With a variable focus wide-angle 25mm UV lens, investigators can scan a large area in a short amount of time. The wide-field, focusable eyepiece assures the elimination of eyestrain.
- Its light-weight body is constructed from a single piece of extruded aluminum which features a User-Facing Control Panel.
- Laser pointer included



Please note: to photograph close-up fine ridge details, you will need the KRIMESITE™ Imager.



No. KSS8900 Scan-n-Find™ Pocket Imager Kit \$9,263.25

Cyclops Imaging System

Specifically designed with crime scene search in mind, the Cyclops RUVIS system unit can be hand-held or configured in a head mount assembly for hands-free use or used in a detachable pistol grip. You'll have comfortable viewing with soft rubber eyecups. The system is fashioned from durable injection-molded material and equipped with a 25mm f/2.8 objective lens and two 15mm easily-focused eyepiece assemblies.



Head Mount



Pistol Grip



Hand-Held



No. KSS7900 Cyclops Imaging System \$10,298.25

RUVIS Systems

KRIMESITE™ Imager Master RUVIS Kit with Cyclops

Designed for agencies that need to work in the field and in the lab, this kit includes both a KRIMESITE™ Imager and a Cyclops unit. It allows investigators to work efficiently alone or as a team to use the Cyclops for scanning the full scene to identify areas of interest and the KRIMESITE™ to capture images of the print with a mounted digital camera (not included). The perfect combination for crime scene technicians, latent print examiners, and for crime laboratories!



Obtain the greatest efficiency at crime scene processing using the KRIMESITE™ and Cyclops Imagers. (Both are contained in the Master RUVIS Kit.)



No. KSS60MRK Master RUVIS Kit with Cyclops \$22,200.75

KRIMESITE™ Master RUVIS Kit with Scan-n-Find

This kit is an alternative to the Master RUVIS kit where the Scan-n-Find takes the place of the Cyclops unit. Some investigators prefer the single eye form and function of the Scan-n-Find™. The kit offers the same efficiency as the other Master RUVIS kit in a different package. This kit features:

- The KRIMESITE™ Imager can be used in the laboratory and for photographing untreated latent fingerprints
- The Scan-n-Find™ is ideal for the crime scene and helps find latent print evidence that is easily overlooked.




No. KSS60SNF Master RUVIS Kit with Scan-n-Find \$22,718.25

RUVIS Systems

Krimesite Capture System

The most versatile RUVIS kit made today. Use in the field and in the lab, with an incredibly easy to use software package and Capture platform to make your search and photo documentation mission simpler.

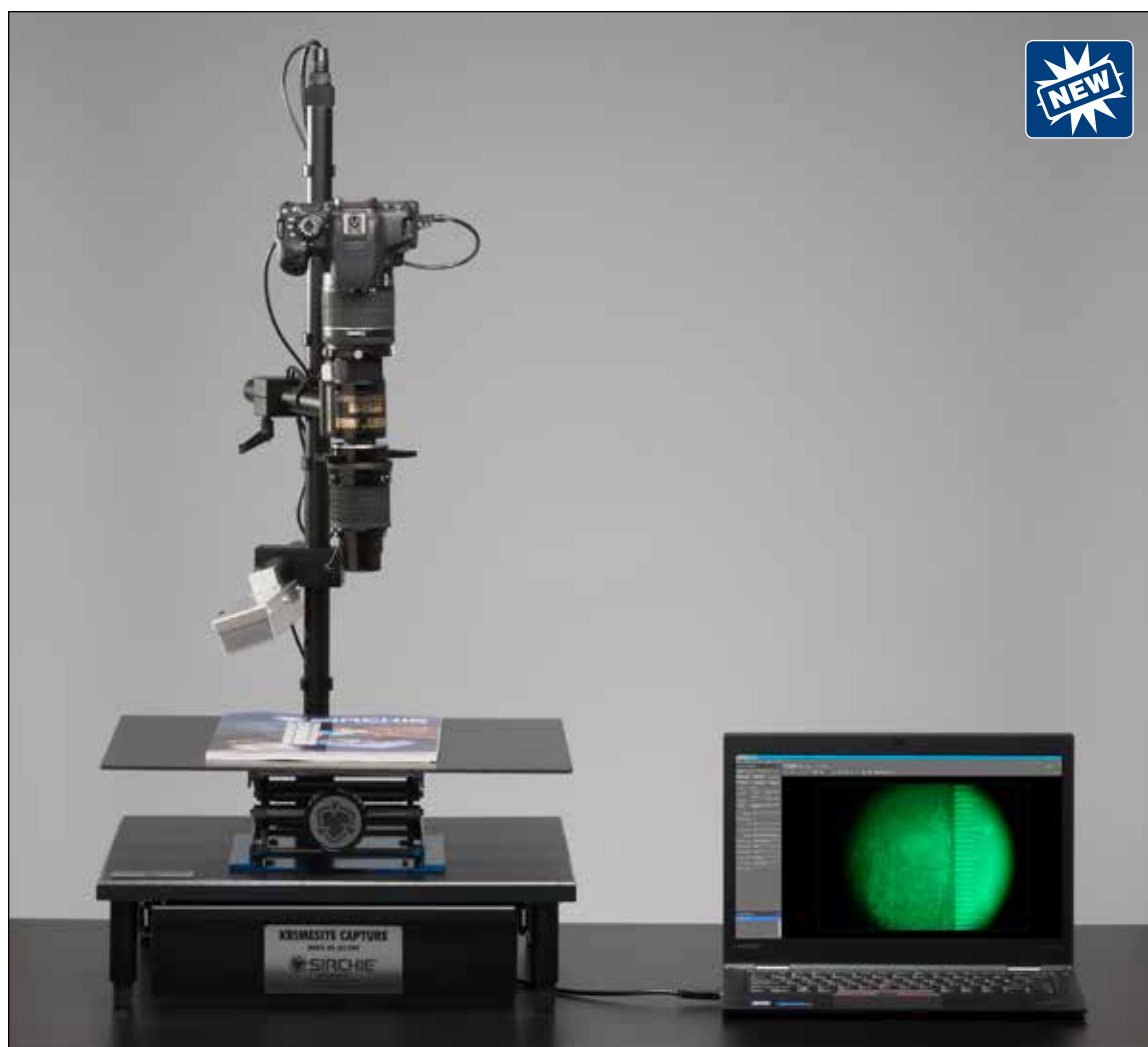

Features include:

- The specialized stand holds the Krimesite Imager and a camera while you use the included software and laptop to capture photographs. You'll be able to control the settings and focus of the camera through the laptop, making your job easier and more efficient.
- User-friendly software allows you to begin looking for latent prints minutes after set-up.
- An adjustable platform gives you the flexibility to photograph evidence on objects of any shape or a range of sizes. Additionally, the imager and camera are fixed on a rotating arm of the stand so you can take photos of large objects.
- A UV light mounted on a flexible arm lets you operate hands-free so you can focus on the image capture.
- Images can be captured in RAW and JPG formats (*RAW is recommended for evidence photography*).

Master Krimesite Capture System available with Nikon or Canon camera

Already have a Krimesite or other RUVIS unit? The Krimesite Capture Upgrade Kit is perfect for you. Also includes Nikon or Canon camera.

Chemical Fingerprint Development

Chemicals, such as cyanoacrylate or ninhydrin, can develop latent prints on different surfaces, develop older prints, or enhance prints already developed by powder so you can see more ridge detail. Using these chemicals can give you significant advantage when developing and preserving latent print evidence.

Latent Print Standards Pad

The Latent Print Standards Pad is an excellent resource to verify that chemicals being used for development of latent prints are working properly so that you can be confident that you have the best chance of developing reliable prints.

This pad is designed to be used for all chemical tests that involve sebaceous oils, amino acids, or salts including:

- Iodine
- DFO
- Ninhydrin
- Silver Nitrate
- Physical Developer



Iodine result prints using the LPSP200 Latent Print Standards Pad



Ninhydrin result prints using the LPSP200 Latent Print Standards Pad



Silver Nitrate result prints using the LPSP200 Latent Print Standards Pad

No. LPSP200

Latent Print Standards Pad

\$11.75

Evidence Collection Course



Our popular Evidence Collection Accelerated Training Class gives you hands-on experience in limited class sizes covering the scientific methods of collection, identification, evaluation, and preservation of physical evidence. Our students love the teach-do method of instruction and take their knowledge back to the field as a lifelong Sirchie graduate.




For more information about this and all other training courses, visit www.sirchie.com or call us at (800) 356-7311.

Chemical Fingerprint Development

Cyanoacrylate Fuming

Super Glue Vapor, ethyl or methyl cyanoacrylate chemically reacts with most latent fingerprints to produce a white, hard deposit of the print.

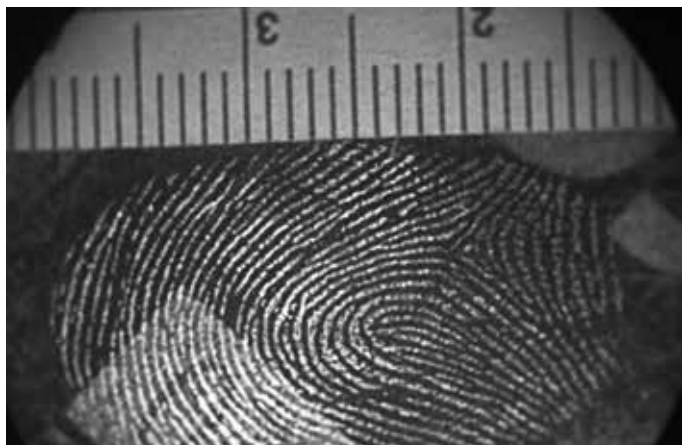

CyanoPowder

Unlike cyanoacrylate liquid which can harden over time rendering it unusable, CyanoPowder has an unlimited shelf life, and can be stored indefinitely. It's easier to handle and safer than liquid glue, designed specifically for lab safety and development performance. Additionally, it's an environmentally friendly alternative to liquid glues, helping labs cut waste and save money. Available in six packs or 30 grams in bulk packs. **Note:** CyanoPowder must be developed at 500° F—use with No. HP110 Fuming Hotplate or No. FR200DT Dual Temp Fuming Chamber.



Often Bought Together:
(CNA106) Disposable Fuming Trays,
(HP110) Hotplate, (FR200DT) Dual
Temp Fuming Chamber

No. CYP15 **CyanoPowder, 1.5g singles, 6 ea.** **\$20.95**

Omega-Print Cyanoacrylate Fuming Compound

Ideal for smaller departments that don't need a large quantity of cyanoacrylate, this compound is supplied with a drop dispensing cap that is resealable to ensure the integrity of the solution while in storage.



No. CNA102 **Cyanoacrylate Fuming Compound, 20g** **\$5.50**

Omega-Print Fingerprint Developer

The leading choice of many law enforcement agencies throughout the world, OMEGA-PRINT comes with four applicator bottles to maintain the integrity of the bulk container. Ideal for departments that process a large number of latent prints.



No. CNA103 **Cyanoacrylate w/4 ea. applicator bottles, 16 oz.** **\$54.95**

Omega-Print Dispersal Pads

Use these pads with any cyanoacrylate liquid fuming compounds.



Often Bought Together:
(CNA103) Omega-Print
Fingerprint Developer



No. CNA104 **Omega-Print Dispersal Pads, 100 ea.** **\$11.25**

Omega-Print Disposable Fuming Trays

Trays are designed to accelerate fuming with cyanoacrylate. Use these trays with or without dispersal pads.



Often Bought Together:
(CNA103) Omega-Print
Fingerprint Developer



No. CNA106 **Disposable Fuming Trays, 25 ea.** **\$4.45**

No. CNA106C **Disposable Fuming Trays, 100 ea.** **\$10.88**

Chemical Fingerprint Development

Iodine Fuming Products

The Fumette Disposable Iodine Fuming Gun

Iodine fuming is very effective for fuming walls and woodwork. The Fumette's design lets you fume vertical and horizontal surfaces and is excellent for small objects. You can control the amount of iodine vapor the whole time. Once you're finished, simply dispose of the whole assembly.



No. DF2016 Fumette Iodine Fuming Gun, 6 ea. \$36.95

Search Iodine Crystal Ampoules

These ampoules are perfect for storing in kits because there's no vapor loss possible. Sealed in glass ampoules, simply break at pre-scored neck when you want to use.



No. A211C Iodine Crystal Ampoules, 6 ea. \$29.95

Ninhydrin Products

Note: Use ninhydrin after iodine and before silver nitrate

The ninhydrin process is ideal for latent prints and blood-contaminated prints on porous surfaces. It's simple to use and very effective. It'll detect prints months or years old.

Ninhydrin Spray

The most popular ninhydrin formulation in the world, it has developed prints over two-years-old in actual cases and 50-year-old prints in experimental cases. Available in 6 oz and 16 oz aerosol cans.



Often Bought Together:
(DFC100) Development Control Chamber, (214CA) Fuming Cabinet & Heat Chamber



No. 201C Ninhydrin Spray, 16 oz. aerosol \$19.95

No. 202C Ninhydrin Spray, 6 oz. aerosol \$9.95

Ninhydrin Crystals

Perfect for agencies that mix their own formulations.




No. NRP01A Ninhydrin Crystals, 25g jar \$26.95

No. NRP02B Ninhydrin Crystals, 100g jar \$71.75

No. NRP03C Ninhydrin Crystals, 500g jar \$284.95

DFO (1,8-Diazafluoren-9-one)

DFO is a fluorescence agent that attaches to latent print residues and causes them to emit light when illuminated with a UV light or alternate light sources. It's ideal for using after ninhydrin.

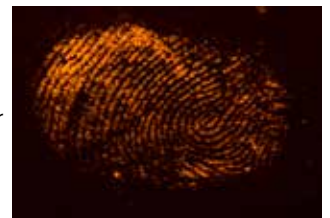

Often Bought Together:
(DFC100) Development Control Chamber, (DFC200) Development Control Chamber



No. DFS300 DFO Pump Spray, 100ml \$39.95

No. LV500 DFO Powder, 1g \$69.95

No. LV5001 DFO Powder, 5g \$315.00

No. LV500L DFO Powder, 50g \$2,795.00

1,2-Indanedione

1,2-Indanedione is a low cost substitute for DFO. Like DFO and ninhydrin, 1,2-Indanedione reacts to amino acid residues on porous surfaces. Developed prints can be difficult to see in normal light, but strongly fluoresce when exposed to 455nm blue light viewed through an orange barrier filter.



No. LV508 1,2-Indanedione, 1g \$29.95

No. LV5081 1,2-Indanedione, 10g \$264.95

Chemical Fingerprint Development

Silver Nitrate

Silver nitrate is especially useful for developing prints on raw wood or when used sequentially on paper documents as the final process after iodine and ninhydrin. It reacts with chloride and produces a silver chloride print.



Search Silver Nitrate Spray

This is a general purpose nitrate solution which is easy to use in the field or the laboratory. Once it's sprayed evenly on the surface, it should dry completely (about two minutes). After drying, use a high-intensity light to view fingerprints.



Often Bought Together:
(SNA4) Silver Nitrate
Accelerator Lamp



No. 205C

Silver Nitrate Spray, 8 oz.

\$29.95

Silver Latent Print Spray

This solution is especially useful when processing documents or objects printed with water soluble inks or dyes.



Often Bought Together:
(SNA4) Silver Nitrate
Accelerator Lamp



No. SLPS300

Silver Latent Print Spray, 8 oz.

\$42.95

Silver Nitrate Crystal Ampoules

These Nitrate Crystals are sealed in glass ampoules for easy dispensing.



Often Bought Together:
(SNA4) Silver Nitrate
Accelerator Lamp




No. SN208

Silver Nitrate Crystal Ampoules, 6 ea.

\$34.95

Silver Nitrate Accelerator Lamp

This shortwave UV lamp rapidly develops latent prints treated with Silver Nitrate and is ideal for both laboratory and field use. **Note:** Use UV protective eyewear when using short wave UV lamps.

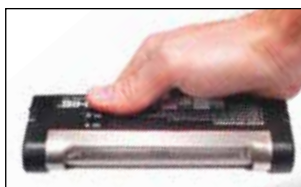

No. SNA4

Silver Nitrate Accelerator Lamp

\$54.95

Physical Developer

Physical Developer Set

Normally, Physical Developer is used on latent prints after they've been treated with ninhydrin. It's actually more sensitive than many Silver Nitrate products.



Available in
1 liter or 3.8 liter sizes.

No. LPD100

Physical Developer Set, 1 liter

\$62.95

No. LPD101

Physical Developer Set, 3.8 liter

\$199.95

Maleic Acid

An acid to use prior to Physical Developer because it neutralizes alkaline binders and fillers. Alkaline will react strongly with the Physical Developer.



Often Bought Together:
(LPD100) Physical
Developer Set

No. LPD110

Maleic Acid, 50g

\$11.75

No. LPD10L

Maleic Acid, 1 liter

\$25.95

Fuming Chambers



It's always better to under-develop prints in a chamber. Overdeveloped prints lose details.

Dual Temperature Cyanoacrylate Fuming Chamber

Our newly redesigned Chamber has a large cubic-foot volume for processing a large number of items simultaneously. The two heating elements can be switched between Standard (230 degrees F) for fuming with regular cyanoacrylate liquid or hi-temp (500 degrees F) to use with our CyanoPowder Crystals. Available in 110V and 220V.



Often Bought Together:
(CNA106) Disposable Fuming Trays,
(CYP15) CyanoPowder, (CNA110)
Cyano-Bloc

No. FR200DT

Dual Temp Cyanoacrylate Chamber, 110V AC

\$1,395.00

No. FR200DT220 Dual Temp Cyanoacrylate Chamber, 220V AC

\$1,450.00

Ductless Drying Cabinets & Accessories

Ductless fume hoods, work stations, DrySafe™ evidence drying cabinets and CyanoSafe™ cyanoacrylate fuming chambers are important tools for any lab. As ductless solutions, these items can be placed anywhere in your lab space to give you flexibility in your layout. They're also popular options in our mobile crime lab vehicles.

DrySafe™

DrySafe™ Evidence Drying Cabinet

We designed this cabinet to protect the operator from harmful bacteria, viruses and putrid odors while ensuring the integrity of your evidence by protecting it from cross-contamination. We also offer a decontamination spray to disinfect the cabinet between cases. Filters are included with all DrySafe™ evidence drying cabinets. See below for available sizes. Available in 110V and 220V.



No. ACEVD30A	DrySafe™, 30" x 28" x 84"	\$6,450.00
No. ACEVD36A	DrySafe™, 36" x 28" x 84"	\$6,850.00
No. ACEVD48A	DrySafe™, 48" x 28" x 48"	\$7,650.00
No. ACEVD60A	DrySafe™, 60" x 28" x 84"	\$10,500.00
No. ACEVD72A	DrySafe™, 72" x 28" x 84"	\$11,800.00
No. ACEVD96A	DrySafe™, 96" x 28" x 84"	\$12,850.00

Visit sirchie.com for full dimensions

DrySafe™ Start-up Kit

We put together this kit so you have everything you need to start using the DrySafe™ evidence drying cabinet.



No. ACEVD34	DrySafe™ Startup Kit	\$245.00
No. ACEVDCLSEP	Clothing Separator, 23" W x 58" L x 1/8" thick	\$175.00

Replacement Filters

Bonded Carbon-Activated Filter

Use this filter to absorb the fumes and odors associated with drying bloody articles.



HEPA Filter (ACFHEPA)

Highly efficient filter for particles, such as toxic powders, asbestos, lead, and fingerprint powders, etc.

Particulate Pre-Filter (ACFPRE)

This filter removes sub-micron particulates to filter out contamination.

No. ACF300	Bonded Carbon-Activated Filter for DrySafe	\$550.00
No. ACFHEPA	HEPA Filter for DrySafe	\$275.00
No. ACFPRE	Pre-Filters for DrySafe, 12 ea.	\$165.00

Kraft Paper Liners

Federal guidelines recommend that you should place a new paper liner on the cabinet bottom before each use. You should retain this liner for potential trace evidence.



No. ACEVD351	Kraft Paper Liners, 100 ea.	\$36.06
--------------	-----------------------------	---------

Kraft Evidence Bags and Paper Roll

Available in four sizes; for larger items, use the paper roll.



No. EB000P	Kraft Paper Bags, 5" x 9 7/8" x 3 1/8", 100 ea.	\$22.34
No. EB001P	Kraft Paper Bags, 7" x 13 3/4" x 4 1/2", 100 ea.	\$24.47
No. EB002P	Kraft Paper Bags, 8" x 18" x 5", 100 ea.	\$26.60
No. EB003P	Kraft Paper Bags, 12" x 18" x 7", 100 ea.	\$31.03
No. KRP100	Kraft Paper Roll, 36" x 1200'	\$87.36

Drying Cabinet Accessories & Fuming Chambers

SearchBreathe Breathable Evidence Bags


This pouch is ideal for evidence items that need breathing capability and visibility. It's ok if the evidence isn't dry when you store it in the pouch.



No. BEP812P	Breatheable Evidence Pouch, 8" x 12", 100 ea.	\$116.20
No. BEP1015P	Breatheable Evidence Pouch, 10" x 15", 100 ea.	\$190.82

Evidence Tags

Federal guidelines recommend using colored tags to indicate drying cabinet status. Blue represents a clean cabinet and Red represents cabinet in use. Space is provided for writing information. Length is 12"

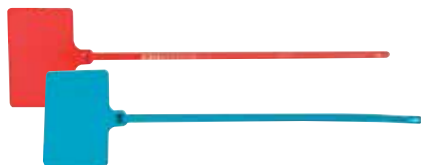

No. ACEVDRT	Red Evidence Tags, 25 ea.	\$39.34
No. ACEVDBT	Blue Evidence Tags, 25 ea.	\$39.34

Solucide Disinfectant Solution

Solucide is an EPA-registered high-level disinfectant and deodorizer spray recommended for cleaning the evidence drying cabinet between cases. FDA approved.



Often Bought with:
(SF0078) Nitrile Gloves

No. ACEVD171	Solucide, 16 oz, 1 ea.	\$16.75
No. ACEVD17	Solucide, 16 oz, 12 ea.	\$190.00

CyanoSafe™

CyanoSafe™ Fuming Chambers and Mini-CyanoSafe™

These fuming chambers are designed to safely process items using the cyanoacrylate fuming method. It provides a controlled environment for processing latent fingerprints on nonporous surfaces while eliminating your exposure to hazardous cyanoacrylate fumes. It features AirSafe, an automatic safety controller of the filter conditions. Available in a benchtop or free-standing design, in 110V or 220V.



Filters are included. Replacement filters available. Refer to owner's manual for replacement schedule.



No. CAS24	mini-CyanoSafe™, Benchtop, 24"	\$4,195.00
No. ACF10018	mini-CyanoSafe™ Hybrid Bonded Carbon Filter, ea.	\$325.00
No. CAS30	CyanoSafe™, Benchtop, 30"	\$5,900.00
No. CAS30T	CyanoSafe™, Free-Standing, 30"	\$8,335.66
No. CAS48	CyanoSafe™, Benchtop, 48"	\$6,750.00
No. CAS48T	CyanoSafe™, Free-Standing, 48"	\$9,500.00
No. CAS60T	CyanoSafe™, Free-Standing, 60"	\$13,500.00
ACF100H	Hybrid Bonded Carbon Filter for CyanoSafe™, ea.	\$525.00
ACFPRE	Pre-Filters for CyanoSafe™, 12 ea.	\$165.00

Visit sirchie.com for full dimensions

CyanoSafe™ Start-Up Kits

The ideal kit to get everything you need to start using your CyanoSafe Fuming Chamber. Available for CyanoSafe and Mini-CyanoSafe.



No. CAS10	CyanoSafe™ Start-Up Kit	\$259.95
No. CAS11	mini-CyanoSafe™ Start-Up Kit	\$89.95

Ductless Forensic Workstations

Dusting Stations

Down Draft Dusting Station

This station gives you or your team critical protection from potentially harmful powders and particulates associated with latent print processing. The station pulls air away from you in two directions: downward and to the rear of the workstation, keeping print powder particulates away from you. It's also easy to maintain and clean and it's light enough to be moved anywhere (around the lab or to crime scenes).



Filters are included. Replacement filters available. Refer to owner's manual for replacement schedule.



No. DDS24	Down Draft Dusting Station, 24" x 17 1/2" x 24"	\$2,400.00
No. DDS36	Down Draft Dusting Station, 36" x 17 1/2" x 24"	\$2,950.00
No. DDS48	Down Draft Dusting Station, 48" x 17 1/2" x 24"	\$3,400.00
No. ACFHEPA8	24" Replacement HEPA Filter	\$158.00
No. ACFHEPA18	36" Replacement HEPA Filter	\$245.00
No. ACFHEPA22	48" Replacement HEPA Filter	\$255.71

Latent Print Dusting Station

Our stations provide superior protection to you or your team while dusting items with a variety of latent print powders. Streamlined air flow pulls particulates away from your breathing zone into a filtration zone. This air flow works because:

- Print powders have a variety of particle sizes and the only effective way to capture the smallest to the largest particles is by pulling them away from your breathing zone.
- Lighter and smaller particles rise while heavier ones will fall to the base and can be cleaned.
- To provide personal protection to you, our station uses an upward airflow
- Our dusting stations have been independently tested for powder containment.



32" Workstations

No. DS100A	Latent Print Dusting Workstation, 32" x 42" x 24"	\$2,750.00
No. ACALGFL	Optional External Fluorescent Light Kit*	\$135.00
No. ACA1011	Optional Steel Cart, 34" fixed height	\$750.00
No. ACFHEPA	HEPA filter for DS100A, 1 ea.	\$275.00
No. ACFPRE1	Pre-Filters for DS100A, 1 ea.	\$185.00

* Light kit must be ordered at time of purchase.

48" Workstations

No. DS200A	Latent Print Dusting Workstation*, 48" x 42" x 24"	\$3,400.00
No. ACA1039	Optional Steel Cart, 34" fixed height	\$895.00
No. ACFHEPA18	HEPA Filter for DS200A, 1 ea. (2 required)	\$245.00
No. ACFPRE4	Pre-Filters for DS200A, 12 ea. (2 required)	\$185.00

* Light kit is included on 48" Workstations.

Ductless Forensic Workstations

Forensic Workstations

Basic Forensic Workstations

Ideal for smaller departments and labs to protect your team against the dangerous particulates of dust and fumes from fingerprint powders, cyanoacrylate, iodine, Ninhydrin, DFO, and drug identification using solvents. Additionally, it can be fitted with a bonded carbon-activated filter for gases or a HEPA filter for particulate containment. Available in 110v or 220V.

32" Workstations



No. AC632C	Carbon Workstation, 32"	\$3,440.00
No. AC632H	HEPA Workstation, 32"	\$3,330.00
No. ACALGFL	Optional Fluorescent Light for 32" Workstation*	\$135.00
No. ACA1011	Optional Steel Cart, 34" fixed height	\$750.00
No. ACF100	Bonded Carbon Filter, 1 ea.	\$425.00
No. ACFHEPA	HEPA Filter, 1 ea.	\$275.00
No. ACFPRE	Pre-Filters, 12 ea.	\$165.00

* Light kit must be ordered at time of purchase.

48" Workstations

No. AC648C	Carbon Workstation*, 48"	\$4,935.00
No. ACA1039	Optional Steel Cart, 34" fixed height	\$895.00
No. ACF10018	Bonded Carbon Filter, 1 ea. (2 required)	\$325.00
No. ACFHEPA18	HEPA Filter, 1 ea. (2 required)	\$245.00
No. ACFPRE4	Pre-Filters, 12 ea. (2 required)	\$185.00

* Light kit is included on 48" Workstations.

Advanced Forensic Workstations

These workstations are for departments that process a large volume of evidence with powder or chemicals and don't want to change filters between applications. This system uses a dual filter bed and these workstations meet and exceed OSHA and ANSI Z9.5 standards and guidelines. Available in 110V or 220V.



32" Workstations

No. AC6200A	Advanced Workstation, 32"	\$4,325.00
No. ACALGFL	Optional External Fluorescent Light Kit*	\$135.00
No. ACA1011	Optional Steel Cart, 34" fixed height	\$750.00
No. ACF100	Bonded Carbon Filter, 1 ea.	\$425.00
No. ACFHEPA	HEPA Filter, 1 ea.	\$275.00
No. ACFPRE1	Pre-Filters, 12 ea.	\$185.00

* Light kit must be ordered at time of purchase.

48" Workstations

No. AC6400A	Advanced Workstation*, 48"	\$6,135.00
No. ACA1039	Optional Steel Cart, 34" fixed height	\$895.00
No. ACF10018	Bonded Carbon Filter, 1 ea. (2 required)	\$325.00
No. ACFHEPA18	HEPA Filter, 1 ea. (2 required)	\$245.00
No. ACFPRE4	Pre-Filters, 12 ea. (2 required)	\$185.00

* Light kit is included on 48" Workstations.

Fingerprint Magnifiers

Magnifiers

Ridge Pro Plus Digital Magnifier

This powerful, hand-held digital magnifier offers a 5-inch wide-screen LCD screen with up to 14x magnification for easy examination. Features include the ability to capture and send an image via USB port, store up to 80 images, built-in stand & lights, built-in handle for field use, and more.



No. FXM50A	Ridge Pro Plus Digital Magnifier	\$1,034.84
------------	----------------------------------	------------

Professional Fingerprint Magnifier



Our most popular magnifier, the JC100P is sturdy, rugged and built for years of use. An optional locking ring will ensure your setting is maintained. Precision focus with four available classification discs: The Henry Disc, the Battley Disc, the Henry/Battley Disc, and the Miracode Palm Disc.



No. JC100P	Professional Fingerprint Magnifier	\$127.87
No. JC101LR	Locking Ring	\$10.38
No. JC101H	Henry Disc	\$26.51
No. JC101B	Battley Disc	\$57.10
No. JC101HB	Henry/Battley Disc	\$51.09
No. JC101M	Miracode Disc	\$51.63

M1 and M2 Fingerprint Magnifiers

Our M1 and M2 magnifiers have superb image clarity and sharpness and feature fine focusing and a lock ring. Both magnifiers provide 4.5x magnification while the M2 features a height adjustment knob. Also available are compatible classification discs: The Henry Disc, the Battley Disc, the Henry/Battley Disc, and the Miracode Palm Disc.



M1 Magnifier

M2 Magnifier

No. PFP100	M1 Fingerprint Magnifier	\$84.76
No. PFP200	M2 Fingerprint Magnifier	\$70.31
No. PFH10	Henry Disc	\$25.40
No. PFB11	Battley Disc	\$47.54
No. PFHB12	Henry/Battley Disc	\$49.44
No. PFM13	Miracode Disc	\$50.54

Rack & Pinion Focus Fingerprint Magnifier

Featuring a two-element 4.5X magnification system with a smooth rack-and-pinion focusing mechanism, this system will ensure you have sharp clear images. Also available are compatible classification discs: The Henry Disc, the Battley Disc, the Henry/Battley Disc, and the Miracode Palm Disc.



No. JC200	Rack & Pinion Focus Magnifier	\$191.84
No. JC201H	Henry Disc	\$39.29
No. JC201B	Battley Disc	\$57.10
No. JC201HB	Henry/Battley Disc	\$51.09
No. JC201M	Miracode Disc	\$51.63

Fingerprint Magnifiers

Ridge PRO Fingerprint Magnifier

Our premium product in our magnifier line, the Ridge PRO is lightweight for its size but rugged and built for hard use. Its textured rubber barrel makes it easy for you to maneuver and with a rubber guard, it's comfortable for the eye and face. The magnifier has a 3X magnification and a 2" wide lens that makes it virtually distortion-free.



No. RPM350	Ridge PRO Fingerprint Magnifier	\$634.30
-------------------	----------------------------------------	-----------------

Super Wide-Field Fingerprint Magnifier

Ideal for fingerprint comparison, this 3.5X magnifier features an exceptionally wide field-of-view. Using the convenient thumbscrew locking adjustment, you can focus in on details quickly. Also available are compatible classification discs: The Henry Disc, the Battley Disc, the Henry/Battley Disc, and the Miracode Palm Disc.



No. JC400	Super Wide-Field Fingerprint Magnifier	\$149.20
No. JC401H	Henry Disc	\$59.55
No. JC401B	Battley Disc	\$72.40
No. JC401HB	Henry/Battley Disc	\$68.79
No. JC401M	Miracode Disc	\$70.47

Aspheric Stand Magnifiers

Not all evidence is on a flat surface so our No. JC5214 4X magnifier features a swivel head for those hard-to-reach areas. Also available in non-swivel head versions in 6X and 8X magnification. All feature a wide undistorted viewing field for less eyestrain, and are lightweight and ready for years of use – also perfect for classrooms!



JC5214 Features a swivel head



Non-swivel head versions available in 6X and 8X

No. JC5214	4X Aspheric Stand Magnifier	\$120.14
No. JC4206	6X Aspheric Stand Magnifier	\$79.71
No. JC4208	8X Aspheric Stand Magnifier	\$90.63

Ridge Counters

Reversible Ridge Counter

A long-time favorite is back by popular demand! This indestructible ridge counter features a reversible point that stows neatly inside the body when not in use.



No. RC221	Reversible Ridge Counter	\$7.59
------------------	---------------------------------	---------------

Retractable Ridge Counter

This counter features a push-button retractable point, which can be replaced.



No. GL004	Retractable Ridge Counter	\$15.02
No. GL0051	Replacement Point	\$4.92

Dual Magnification Fixed Focus Magnifier

This magnifier gives you two options of magnification: 3.5 and 7. It's intended to verify fingerprints and documents with security features.



No. FM1005	Dual Magnification Fixed Focus Magnifier	\$83.95
-------------------	-------------------------------------------------	----------------

Forensic Examination Station

Digital Forensic Evidence Examination Station

The Digital Forensic Evidence Examination Station makes evidence examination more efficient and it's easy to use for forensic physical examination, which takes many forms, shapes, and sizes. New 1080HD Display!



You need this station if you examine:

- Tools and toolmarks
- Fingerprint/palm print record cards and latent prints
- Footprint evidence
- Trace evidence: hair, fibers, glass fragments, soil, etc.
- Documents: printed and handwritten
- And much more.

With this station, forensic examiners will have:

- The ability to view small details with a moveable viewing table and fully adjustable magnification – up to 50X.
- A wide range of color correction, brightness, and contrast on an oversized 22-inch LCD, high-definition viewing screen.
- Plenty of room on the mounting table to exam 3D objects. Magnets can also secure evidence to the evidence platform.
- Connect and control from any computer with the FXE350A software, which is included. You can use your own computer to turn the station's lights on and off, adjust screen colors, pause the camera, and adjust the brightness of the images.

No. FXE350A

Digital Forensic Exam Station

\$4,341.83




Fully adjustable magnification (up to 50X) allows you to examine minute details of forensic evidence.





Moveable work table is fully adjustable



User-friendly capture software included

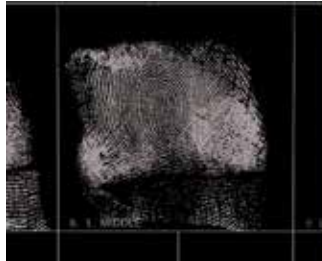

Use onboard controls to adjust brightness, contrast & reverse or colorize picture

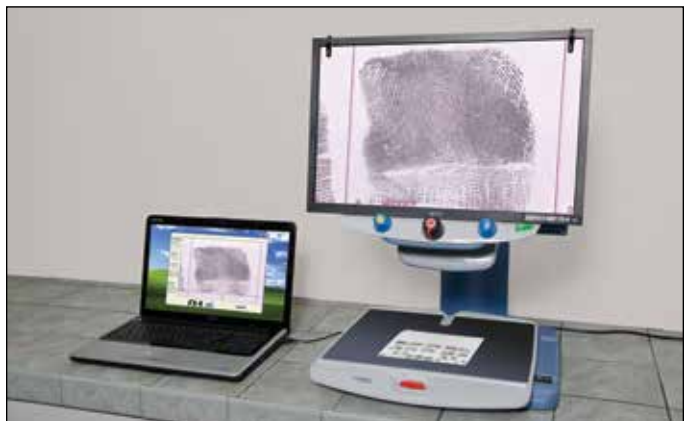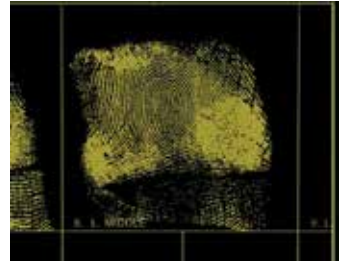

FXE350A Digital Forensic Evidence Examination Station includes Digital Screen Capture. Laptop computer not included.


DNA & Biological Evidence Collection

The analysis of DNA is similar to fingerprints in trying to determine a match. With identifying features, DNA can either link or eliminate a suspect from a crime as well as link crime scenes together. The most frequent forms of DNA found at a crime scene are: blood, semen, hair, saliva, skin cells, mucous, perspiration, urine, fingernails, etc. The more we know how to use and collect DNA, the more powerful a tool it becomes.

Sexual Assault Kits

Sexual Assault Victim Kit

Recently updated with input from SANE-A nurses, this kit is for adult and adolescent examinations and contains all the necessary forms and medical tools to conduct a thorough yet compassionate exam.

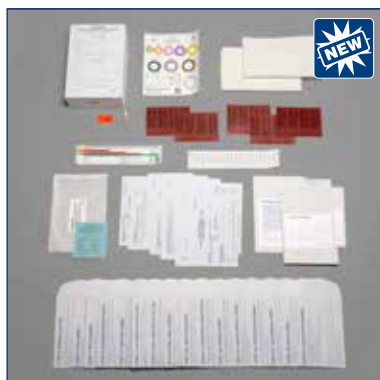


No. VEC200A

Sexual Assault Victim Kit

\$19.95

Pediatric Sexual Assault Kit

Developed with input from SANE-P nurses and the state of Alaska, this sexual assault kit was developed to address the unique needs presented in the examination of a possible sexual assault victim who has not yet entered puberty.

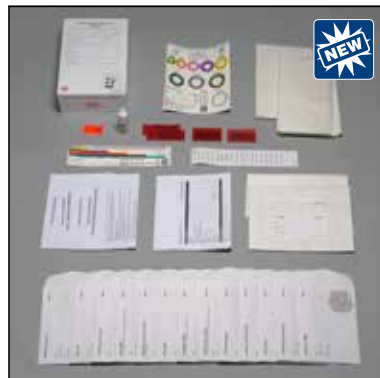


No. PEDSA100

Pediatric Sexual Assault Kit

\$17.08

Suspect Sexual Assault Kit

Proper evidence collection for a possible suspect in sexual assault exams can be key to establish culpability or innocence by linking the perpetrator and the victim through DNA. All necessary components, with clear labeling and step-by-step organization helps simplify the collection procedure.

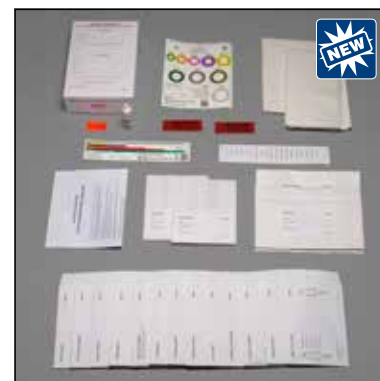

No. SEC200

Suspect Sexual Assault Kit

\$15.01

Custom Sexual Assault and Biological Evidence Collection Kits

As a manufacturer, we have full capabilities to create custom kits for the proper collection of DNA and other biological evidence. We work with jurisdictions and organizations around the world to produce kits built to their exact specifications and provide expertise and guidance on best practices to ensure the highest quality up-to-date products.



Contact us today to discuss your needs and receive a tailored end to end solution delivered with world-class customer service.

DNA & Biological Evidence Collection

DNA Collection Kits

Buccal Swab Kit

With this kit, non-medical personnel can collect buccal cell samples for DNA testing by swabbing the inside of the subject's cheek.



No. BSK100

Buccal Swab Kit

\$3.15

Crime Scene DNA Evidence Collection Kit

You should always be prepared to collect potential sources of DNA and this kit has everything you need in one place. The most frequent forms of this valuable evidence are: fresh blood, semen, hair roots, or clothing stained by physiological fluid.



No. DNA100

Crime Scene DNA Evidence Collection Kit

\$21.80

Suspect Evidence Collection Kit (Buccal Swab/Finger Stick)

Collect evidentiary samples from a sexual assault suspect. Step-by-step instructions guide through collecting and preserving evidence from clothing, hair samples, buccal swabs, and blood.



No. SEC100BF Suspect Collection Kit (Buccal Swab/Finger Stick) \$9.86

Blood & Urine Collection Kits

Blood Specimen Collection Kit

Designed for medical personnel, this kit has everything to collect a whole blood sample for DNA testing, preserve it, and send to the laboratory.



No. BSC50

Blood & Urine Collection Kit

\$5.40

Blood/Urine Specimen Collection Kit

This kit has everything you need for taking blood samples and urine samples and sending to the laboratory. Step-by-step instructions will lead you through documentation to packaging the specimens correctly for shipment.



No. BUK200

Blood/Urine Specimen Collection Kit

\$8.48

Blood Alcohol Collection Kit

Our comprehensive blood alcohol collection kit is designed to provide the safe and efficient collection, handling, and mailing of blood specimens to the lab.



No. BSC100A

Blood Alcohol Collection Kit

\$9.22

DNA & Biological Evidence Collection

Presumptive Seminal Fluid Identification

The Acid Phosphatase Test

In sexual assault cases, you can save time with a presumptive non-destructive chemical test to identify possible seminal stains.

The Acid Phosphatase Test is one of those presumptive tests and it's considered to be one of the most conclusive, non-spermatozoa test for the presence of semen. The reagents for this test are sealed in glass ampoules to ensure integrity and long life. You should conduct a preliminary exam with longwave UV lights or an alternate light source.



TIP: The reagents should not be used directly on the suspect stain. Using the contact method, press a filter paper wetted with distilled water against the stain and then test the paper. This preserves the stain for later confirmatory and possible DNA testing.

Seminal Fluid Dischaps

Easy-to-use with premeasured reagents individually sealed in glass ampoules with polyethylene applicators. This helps ensure against contamination.



No. DCS100 Seminal Fluid Dischaps, 6 ea. \$31.93

Master Seminal Fluid Test Kit for Detection of Acid Phosphatase

A complete master field kit, you'll be able to perform presumptive tests of suspected seminal tests.



No. SF298 Master Seminal Fluid Test Kit \$142.00

RSID Detection Field Test Kits

Easy-to-use tests in the field as confirmatory tests for human fluids. There's no cross reactivity with other bodily fluids or animal fluids. The RSID Bodily Fluid Detection Field Test Kit includes an assortment (10 tests total) of presumptive blood tests, semen tests and saliva tests.



No. IDSE200 RSID Semen Detection Field Test Kit, 10 tests \$217.99

No. IDSA300 RSID Saliva Detection Field Test Kit, 10 tests \$196.14

No. IDBF400 RSID Bodily Fluid Detection Field Test Kit, 10 tests \$207.07

Fingernail Collection Kits

Fingernail/Fingertip Swabbing & Clippings Kit

Collect fingernail swabbings and clippings with this kit; it contains necessary items for each hand as well as evidence seals to ensure evidence integrity.



Fingernail/Fingertip Swabbing & Clippings Kit



Fingernail/Fingertip Swabbing Kit

No. FSWCL100PT Swabbing & Clipping Kit (w/cotton & poly swabs) \$5.19

No. FSW100 Swabbings Kit (w/cotton swabs) \$3.60

No. FSW100PT Swabbings Kit (w/cotton & poly swabs) \$3.81

DNA & Biological Evidence Collection

Fingernail Scrapings Kit

Get evidence from fingernails by using this kit. You can keep evidence separate from each hand with the included swabs, envelopes, and evidence seals.



No. FSC100 **Fingernail Scrapings Kit** **\$2.50**

Fingernail Clippings Kit

Get evidence from fingernail clippings by using this kit. Keep evidence separate from each hand with the included swabs, envelopes, and evidence seals.



No. FCL100 **Fingernail Clippings Kit** **\$3.48**

EZ-Dry™ Swab Boxes and Evidence/Integrity Seals and Kit

The EZ-Dry Swab Box is one of the best ways to collect, package, label, and seal cotton swabs without the possibility of sample switching or secondary contamination. It also provides a way to maintain a continuous chain of possession. Use our new FSD100 forensic dryer to dry swabs.



No. KCP248C **EZ-Dry Swab Boxes, 100 ea.** **\$53.25**

No. SB102C **EZ-Dry Swab Box Evidence/Integrity Seals, 100 ea.** **\$15.94**

No. SB1000 **EZ-Dry Swab Box Kit** **\$7.05**

Fab-Swab DNA Collection Applicator

This convenient, reclosable, dry transport device is designed for reliable sample and evidence collection. The patented tube has a breathable filter label that covers the swab to allow drying while protecting against contamination.



No. FSWAB50 **Fab-Swab DNA Collection Applicator, 50 ea.** **\$65.51**

Cap-Shure® Sterile Forensic Evidence Collection Swabs

These swabs are a unique crafted method for collecting small quantities of very fragile evidence. Whether it's blood, saliva, latent print residue or other physiological fluids, these Cap-Shure® Swabs ensure the integrity of the samples



No. CP100 **Cap-Shure® Sterile Swabs, 100 ea.** **\$102.88**

No. CP500 **Cap-Shure® Sterile Swabs, 500 ea.** **\$479.72**

No. KCP254C **Vented Swab Boxes, 3/4" x 3/4" x 6 1/4", 100 ea.** **\$49.58**

HydraFlock® Sterile Swabs

HydraFlock swabs produce superior sample collection, rapid elution, and increased test sensitivity by combining a unique microstructure with innovative fibers.



No. HFS50 **HydraFlock Sterile Swabs, 50 ea.** **\$20.70**

Purflock® Sterile Swabs

Made with nylon fibers, these swabs have the ability to better absorb organisms at the collection site and during transport.



No. PFS50 **Purflock Sterile Swabs, 50 ea.** **\$18.52**

Swab Drying Rack

Reduce the risk of cross-contamination while applicators are drying with this swab box. It holds 12 applicators and is preprinted with descriptions to keep track of where applicators came from.



Often Bought with:
(KCP160C) Sterile Cotton
Applicator Swabs

No. SDR100C **Swab Drying Rack, 100 ea.** **\$72.44**

DNA & Biological Evidence Collection

DNA Free Specimen Collection

DNA-Free Swabs

Certified to be DNA-Free, these applicators are made with high-quality USP grade cotton.



No. FDNA100 DNA Free Swabs, 100 ea. \$26.60

DNA-Free Ultra Flocked Swabs


DNA-Free swab with a plastic handle and standard tip.



No. FDNA350 DNA Free Ultra Flocked Swabs, 50 ea. \$79.90

SecurSwab™ 2

These DNA-Free swabs come with barcodes to provide the end-user with control in labeling the sample and maintaining chain-of-custody control.

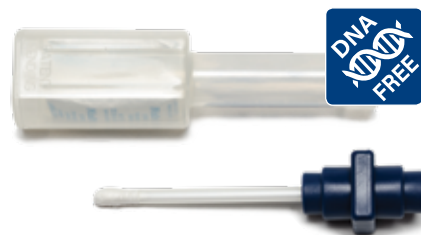

No. SWB100N SecurSwab 2, 100 ea. \$265.45

No. SWB100B SecurSwab 2 w/barcode, 100 ea. \$292.10

DNA Collection Cards

These collection cards are treated with a proprietary chemistry intended to prevent environmentally induced degradation. They allow long term preservation of DNA from biological fluids stored at ambient temperature. Use **non-indicating** cards for blood sampling, and **indicating** cards (that turn color from pink to white) for saliva, buccal cells and urine.



No. DNAG1 DNA Card, 1-Spot (non-indicating), 100 ea. \$239.95

No. DNAG2 DNA Card, 2-Spot (non-indicating), 100 ea. \$305.95

No. DNAG4 DNA Card, 4-Spot (non-indicating), 100 ea. \$375.95

No. DNAG1C DNA Card, 1-Spot (indicating), 100 ea. \$295.95

No. DNAG2C DNA Card, 2-Spot (indicating), 100 ea. \$365.95

No. DNAG4C DNA Card, 4-Spot (indicating), 100 ea. \$435.95

Forensic Swab Dryer

Use our Forensic Swab Dryer to provide a safe, controlled environment for drying swabs after specimen collection. This dryer will hold up to 72 swabs with a drying time of one hour. This unit delivers great drying results in a relatively short period while protecting the integrity of your evidence. Available in 110V and 220V.

No. FSD100 Forensic Swab Dryer, 110V AC \$695.00

No. FSD100220 Forensic Swab Dryer, 220V AC \$695.00



Often Bought with:
(ACEVDRT) Red Evidence Tags,
(SDR100C) Swab Drying Racks



DNA & Biological Evidence Collection

Specimen Collection Tubes and Accessory Products

Sterile & Distilled Water

Pre-moisten applicator swabs when collecting dry specimens with sterile or distilled water. Avoid over-saturating swabs.



Often Bought with:
(KCP160C) Sterile Cotton Applicator Swabs



No. KCP247C Sterile Water, 3ml vial, 100 ea. **\$37.26**

No. DNA1009 Non-Sterile Distilled Water, 1 oz. bottle, 1 ea. **\$2.24**

Blood Collection Tubes

Gray, Lavendar, and Red stopper collection tubes available.



Often Bought with:
(SF0078) Powdered Nitrile Gloves



No. KCP176C Gray Stoppered Blood Tube, 10 ml, 100 ea. **\$76.43**

No. KCP122C Lavendar Stoppered Blood Tube, 10 ml, 100 ea. **\$60.04**

No. CC110C Red Stoppered Blood Tube, 10 ml, 100 ea. **\$58.95**

Liquid Absorption Sheets

Ideal anywhere where bodily fluids are being collected, these sheets are also perfect to ship along with specimen containers during transport.



Often Bought with:
(KCP176C) Gray Stoppered Blood Tube,
(KCP122C) Lavendar Stoppered Blood Tube,
(CC110C) Red Stoppered Blood Tube



No. KCP113C Liquid Absorption Sheets, 2 1/2" x 4", 100 ea. **\$9.78**

Padded Blood Tube Containers

Large enough to hold one 10ml tube.



Often Bought with:
(KCP176C) Gray Stoppered Blood Tube, (KCP122C) Lavendar Stoppered Blood Tube, (CC110C) Red Stoppered Blood Tube



No. BTC100C Blood Tube Container and Plug Insert, 100 ea. **\$68.79**

Blood Tube Containers and Caps

Large enough to hold one 10ml tube and it provides protection for blood tubes.



Often Bought with:
(KCP176C) Gray Stoppered Blood Tube, (KCP122C) Lavendar Stoppered Blood Tube, (CC110C) Red Stoppered Blood Tube



No. BTC200C Blood Tube Containers & Caps, 100 ea. **\$43.66**

Bubble Pouch Tube Protectors

These bubble pouch tube protectors have a peel and seal flap that keeps the tubes secure in the pouch.



Often Bought with:
(KCP176C) Gray Stoppered Blood Tube, (KCP122C) Lavendar Stoppered Blood Tube, (CC110C) Red Stoppered Blood Tube



No. KCP253C Bubble Pouch Tube Protectors, 100 ea. **\$38.18**

Develop every skill.



SIRCHIE[®]
Education & Training

We deliver **hands-on experiences** with innovative and practical exercises to **further your career**. Continue learning through our online webinars and Lunch-n-Learns. Access to this knowledge doesn't end: **contact our experienced instructors** anytime when you need help in the field.



 **SIRCHIE®**
Command Every Scene®

Phone: 800.356.7311 • 919.554.2244

Email: sales@sirchie.com

www.sirchie.com